



ATA1200, ATA2200
ATZ1200, ATZ2200

ATA220÷ATA520
ATZ220÷ATZ520

USER MANUAL

ATA and ATZ Series



Contents:

1.	General description	3
2.	Set	3
3.	Safety rules	4
4.	Technical data	5
5.	General scale view	6
6.	Keys and indicators	8
7.	Preparing working environment	9
8.	Preparing scale to work	10
9.	General operation principles	11
10.	Internal calibration	11
11.	Checking the balance	13
12.	Connection with a computer or a printer	14
12.1	Detailed Long protocol description	15
13.	Start-up	18
14.	Weighing with tare	19
15.	Scale menu	20
16.	Menu navigation rules	21
17.	Scale setup (SEtUP)	27
17.1	External standard calibration and internal calibration options (CALIb)	28
17.2	Autozeroing function (AutotAr)	32
17.3	Weight unit selection (UnIt)	33
17.4	Serial port parameters setting (SErIAL)	34
17.5	Printout configuration (PrInt)	35
17.6	Weighing speed selection (SPEED)	37
18.	Special functions description	38
18.1	Product and user identification (Prod and USEr)	39
18.2	Pieces counting function (PCS)	41
18.3	Weight unit selection (UnIt)	42
18.4	Percentage weighing function (PErC)	43
18.5	Label choosing function (LAbEL)	44
18.6	Weighing animals function (LOC)	45
18.7	Constant tare memory function (tArE)	46
18.8	Maximum value indication function (UP)	48
18.9	Total weight function (totAL)	49
18.10	Checkweighing function (thr)	51
18.11	Setting date and time function (dAtE)	54
18.12	Statistical calculations function (StAt)	55
18.13	Function for summing recipe ingredients (rECIPE)	58
18.14	Density determination (dEnSItY)	59
18.15	Paper weight determination function (PAPEr) - optional	62
19.	Troubleshooting and maintenance	63
	Notes	

1. General description

ATA i ATZ series electronic scales are destined for laboratory works which require high accuracy and for wide range of technical purposes as well. ATA balances have internal calibration, which corrects weighing precision during exploitation. ATZ series are not equipped with internal calibration system.

SPEED option enables to change weighing speed that allows to adjust the balance to working conditions. In particular, if user chooses *FAST* speed setting then small weight portion batching for eg. in pharmacy is possible.

All scales are metrologically tested by manufacturer.

All balances can have legal verification (except ATZ520) or be calibrated by laboratory with PCA accreditation.

According to an order balances can be calibrated.


NACE classification: 33.20.31.

2. Set

Standard set consists of:


1. Scale,
2. Pan and pan support,
3. Feeder,
4. Draft shield with cover (option),
5. User manual,
6. Guarantee card.

3. **Safety rules**

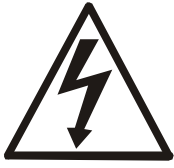

It is necessary to follow safety rules of work with the scale shown below. Obeying those rules is the condition to avoid electrical shock or damage of the scale or connected peripheral devices.

- All repairs and necessary regulations can be made by authorised personnel only.
- To avoid fire risk use a feeder of an appropriate type (if feeder is supplied with the scale) and supply voltage has to be compatible with specified technical data.
- Do not use the scale when its cover is opened.
- Do not use the scale in explosive conditions.
- Do not use the scale in high humidity environment.
- If the scale seems not to operate properly, switch it off and do not use until checked by authorised service.

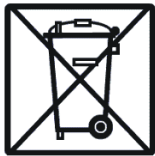

According to current acts of law about protection of natural environment, wasted scales should not be put into waste containers together with ordinary waste.

- Wasted scale after operation period can be delivered to units authorized for gathering wasted electronic devices or to the place where it was bought.

4. Technical data


Type	ATA220 ATZ220	ATA320 ATZ320	ATA520 ATZ520	ATA1200 ATZ1200	ATA2200 ATZ2200
Load (Max)	220g	320g	520g	1200g	2200g
Readout unit (d)	0,001g	0,001g	0,001g	0,01g	0,01g
Verification plot (e)	0,01g	0,01g	0,01g	0,1g	0,1g
Tare range	-220g	-320g	-520g	-1200g	-2200g
Accuracy class	II				
Working temperature	+10 ÷ +40°C				
Weighing time	< 3s			< 2s	
Pan dimension	Ø115mm			Ø150mm	
Dimensions	185x290x90mm				
Interfaces	In standard: RS232C and USB; Options: LAN, Wi-Fi or RS485				
Supply	~230V 50Hz 6VA / =12V 1A				
Scale weight	ATA : 2,6kg ATZ: 2,1kg				
Recommended standard of mass	F2 200g	F2 200g	F1 500g	F2 1000g	F2 2000g

Note:

F2 and F1 are names of international calibration weight classes according to O.I.M.L. Requirements about calibration weight accuracy are connected with these classes.


5. General scale view

ATA120÷ATA520 and ATZ120÷ATZ520 scale:

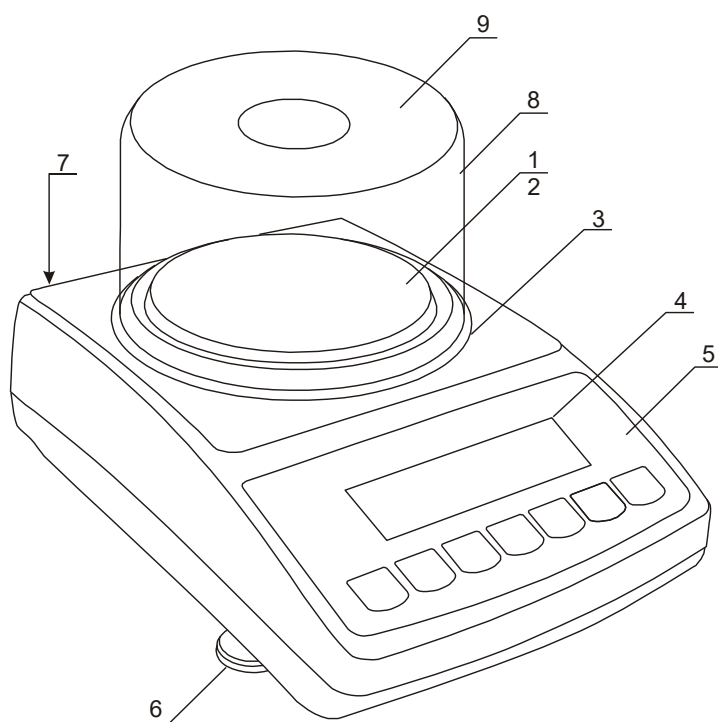


- 1 – pan
- 2 – pan support (under pan)
- 3 – pan ring (against blows)
- 4 – display LCD
- 5 – keys
- 6 – rotating legs
- 7 – water level
- 8 – draft shield (option)
- 9 – draft shield cover (option)

ATA1200÷ATA2200 and ATZ1200÷ATZ2200 scale:

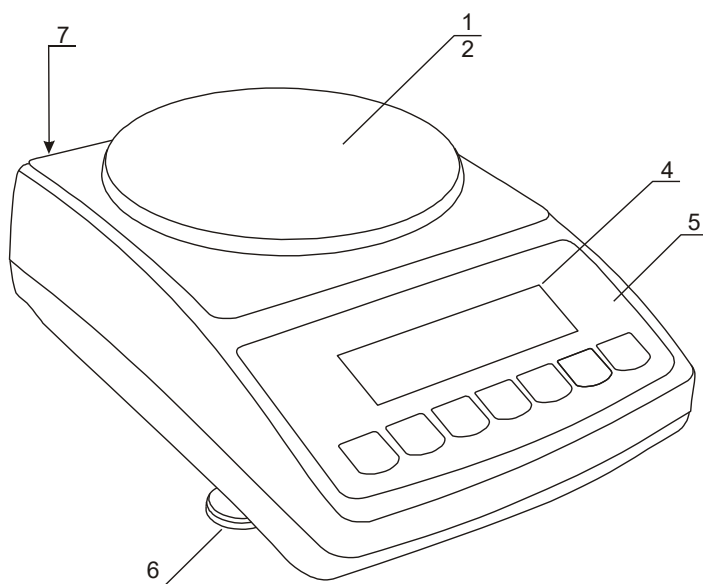


- 1 – pan
- 2 – pan support
- 3 – information window
- 4 – display LCD
- 5 – keys
- 6 – rotating legs
- 7 – level

Connectors view:

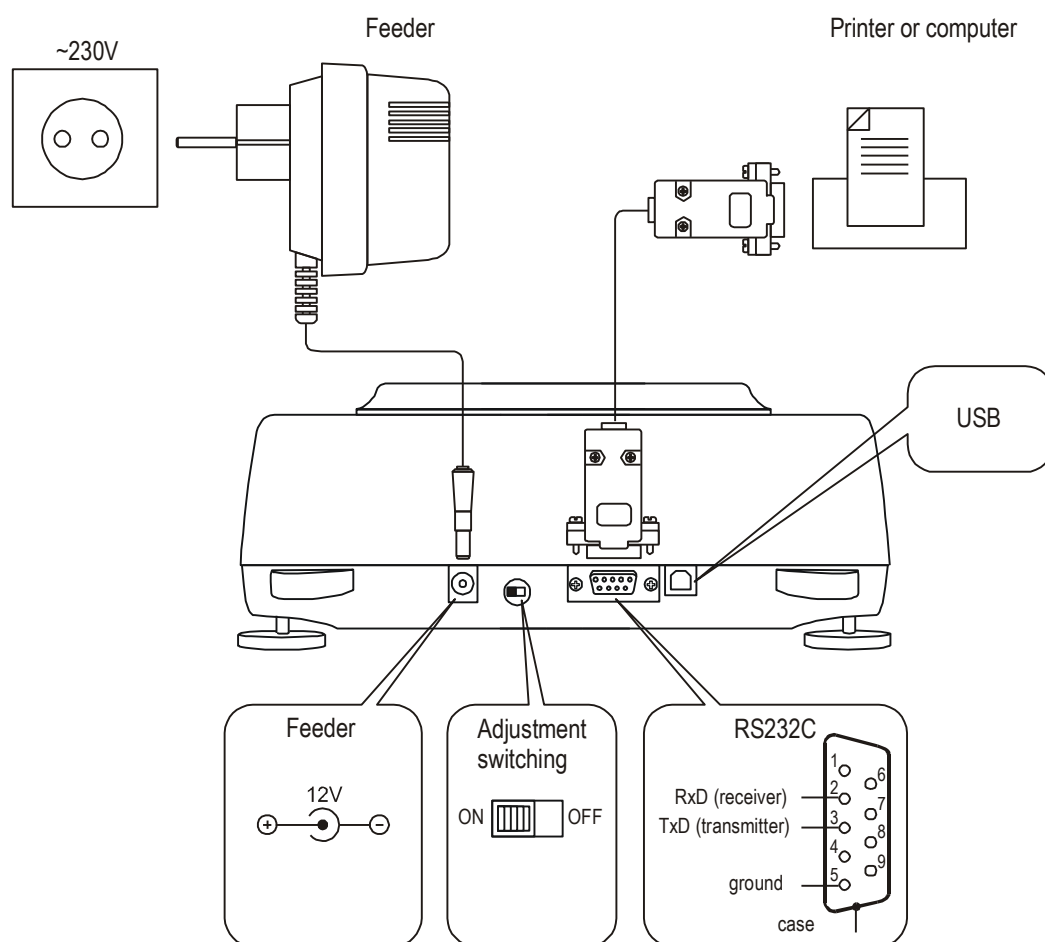

6. Keys and indicators

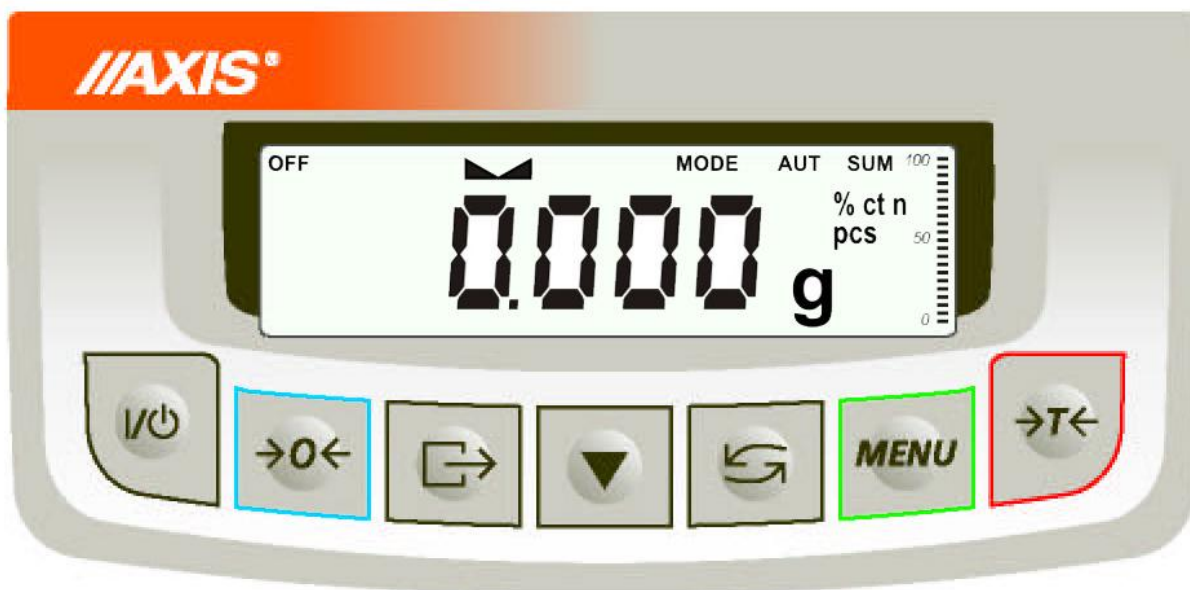

key	I/O	- switch on / switch off (standby),
"	→T←	- taring (storing package mass subtracted from weighed mass)
"	↺↻	- change mode of balance work,
"	→0←	- zeroing the scale when pan is empty (option) ,
"	MENU	-special function menu,
"	□→	- result printout,
"	▼	- internal calibration / changing menu position,
indicator	→0←	- zero indicator (when scale pan is empty),
"	—	- indicator of weighing result stabilisation,
"	NET	- net mass (after use of →T← key),
"	MODE	- indicator of switching special function on,
bar graph indicator		- indicator of scale load (0-100%).
indicator	OFF	- switching scale with I/O key (standby),
"	pcs	- indication in pieces



The use of keys during entering numeric values (special functions):

- ▼ - increment current digit,
- - insert comma,
- T← - move to next position,
- MENU - finish entering

Note:

→0← keys and →0← indicators only work in balances with d=e.


7. *Preparing working environment*

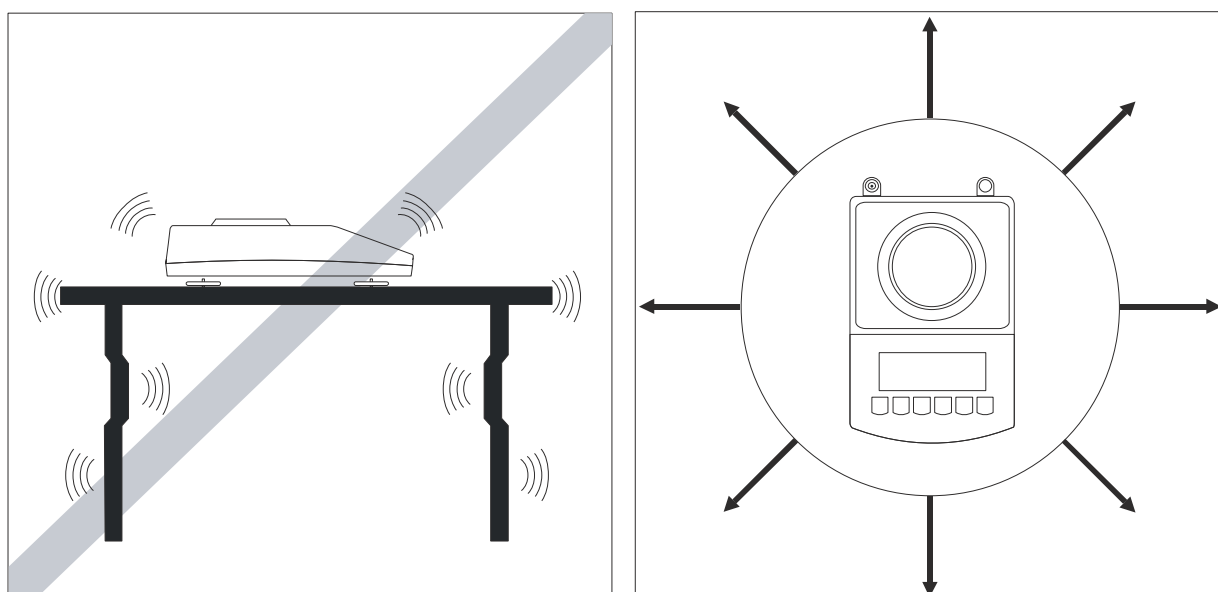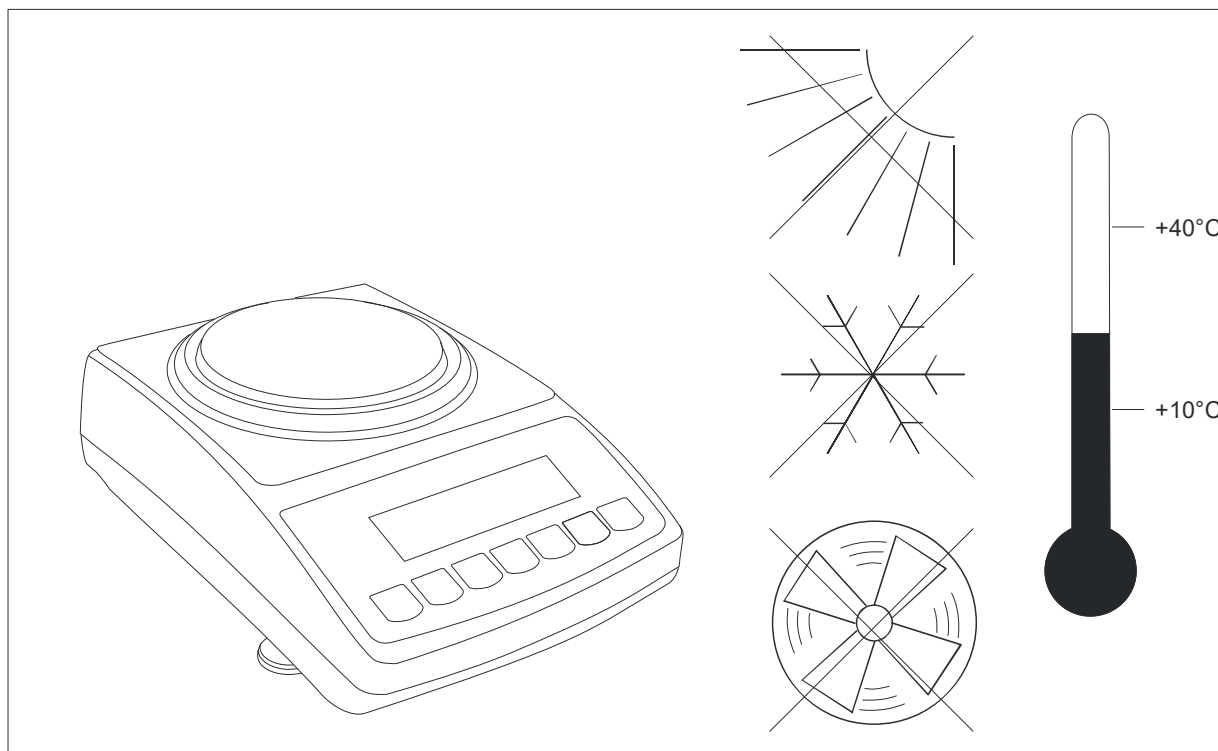

Location for the scale should be chosen with care in order to limit influence of the factors that can interrupt working scale. This location has to maintain proper temperature for working scale and necessary space for its operating. The scale should stay on stable table made of material that does not influence magnetically on the scale.

Rapid air blasts, vibrations, dust, rapid temperature changes or air humidity over 90% are not allowed in scale surrounding. The scale should be far from heat sources and devices emitting strong electromagnetic or magnetic fields.

8. Preparing scale to work

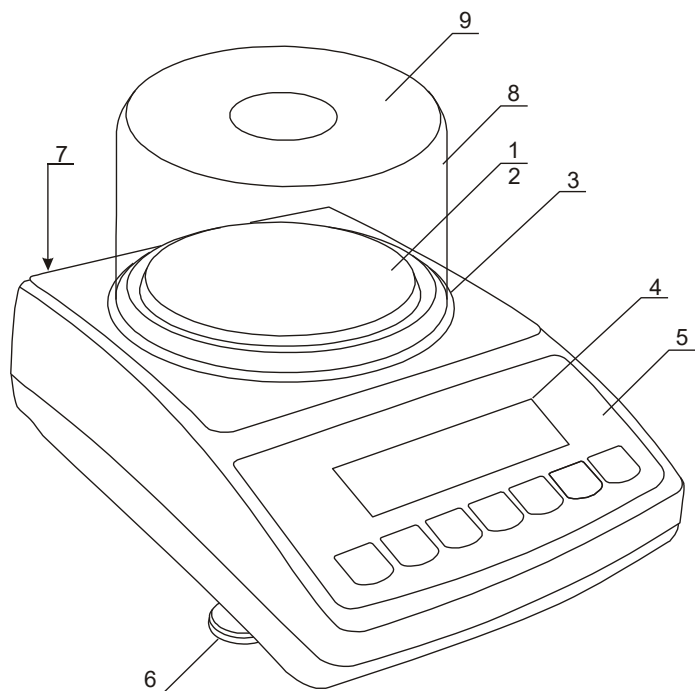

1. Take the scale and feeder out of the package. It is recommended to keep the original scale package in order to transport the balance safely in future.


2. Place the scale on a stable ground not affected by mechanical vibrations and airflows.

3. Level the scale using rotating legs 6 so that the air bubble in water level 7 at the back of the scale is in the middle.


4. Put on the ring securing the weighing pan 3 (applies to scales with a Ø115mm weighing pan).

5. Put draft shield 8 on the scale and cover 9 on it (option).

6. Level the scale using the rotating feet 6 so that the air bubble in the level 7 located at the back of the scale is in the middle position.

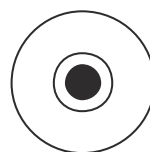

Good

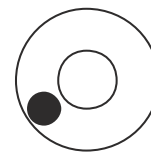

Bad

Leveling a scale equipped with a bubble level involves adjusting the tilt using the rotating feet and observing the gas bubble inside the bubble level. When the bubble is exactly in the center, the scale is level.



If the scale was taken from a lower temperature surrounding to a room with higher temperature, e.g. in winter, moisture can liquefy on the scale casing. Do not connect power supply to the scale, because this can cause damage or improper work of the scale. In this case leave the scale for at least 4 hours unplugged for acclimatization.



After each change of the scale setting, the scale should be levelled and in ATA scales, internal calibration should be performed using the key ▼.

9. General operation principles

1. For ATA scales (equipped with an internal calibration mechanism), it is recommended to perform internal calibration before each series of important measurements during operation or enable automatic calibration.
2. Weighed mass should be placed in the middle of the pan.
3. The scale allows taring in the whole measuring range. To tare the scale press $\rightarrow T \leftarrow$ key. Taring does not extend measuring range, but only subtracts tare value from mass value of a sample placed on the pan. To make the control of a load on the pan easier and to avoid exceeding measurement range, the scales have load indicator calibrated $0 \div 100\%$.
4. Weighing result should be read when the indicator $\blacktriangle \blacktriangle$ lights, which signals result stabilisation.
5. When the scale is not used but it is necessary for it to be ready to work, it can be switched off by pressing I/ϕ key. The scale reading system is then switched off and scale goes to standby mode signalled with *OFF* indicator. Switching the scale on is preformed by pressing I/ϕ key.
6. In sales having $\rightarrow 0 \leftarrow$ key (zeroing) active it should be checked if zero indicator $\rightarrow 0 \leftarrow$ is displayed before sample is placed on the pan. If not, press $\rightarrow 0 \leftarrow$ key and wait until the scale is zeroed and zero indicator appears. After that load can be placed on scale pan.
7. Scale mechanism is a precise device sensitive to overweight, mechanical shocks and strokes.
8. After every change of balance position, level the balance and perform internal calibration.
9. The scale should not be used to weigh ferromagnetic materials due to reduced weighing accuracy.



Do not overload the scale more than 20% of maximum capacity.
Do not press the pan with a hand.



For transportation time, pan support and pan ring should be packed separately.



In order to confirm correctness of the scale during its operation, before starting and after finishing every valid measurement series it is recommended to check weighing accuracy putting calibration weight or other object of exactly known mass on the scale. In the case when allowable measurement error of the scale is exceeded, it is recommended to perform calibration with external weight or contact authorised service centre.

10. Internal calibration

ATA scales are equipped with an internal calibration system (ATZ scales do not have it), the purpose of which is to maintain the required accuracy of the scale under changing operating conditions.



For internal calibration, use the key (press twice). Calibration should be performed before each series of important measurements, as well as after the flashing ▼ CAL icon appears.


To interrupt the internal calibration, press the key during the calibration ▼, which will display the *CAL End* message and allow you to continue the measurements.

Internal calibration involves the automatic application of an internal mass standard by the scale's mechanism and the introduction of an accuracy correction in the scale's software. This correction is necessary due to differences in gravitational acceleration between the scale's location and its place of operation, as well as changes in the scale's level, temperature, etc.

In the first hour after powering on, the scale's internal calibration occurs several times. This is indicated by the *Aut_CAL* message, which appears for one minute before the automatic internal calibration begins (the user can shorten this time by pressing the key twice ▼).

In the following hours, internal calibration occurs after a pre-set operating time and with any ambient temperature change exceeding a pre-set value (for verified scales, this is every 2 hours and 2 °C). The flashing CAL icon indicates when internal calibration is due. Available internal calibration options are provided in section 17.1.

To perform internal calibration, follow these steps:

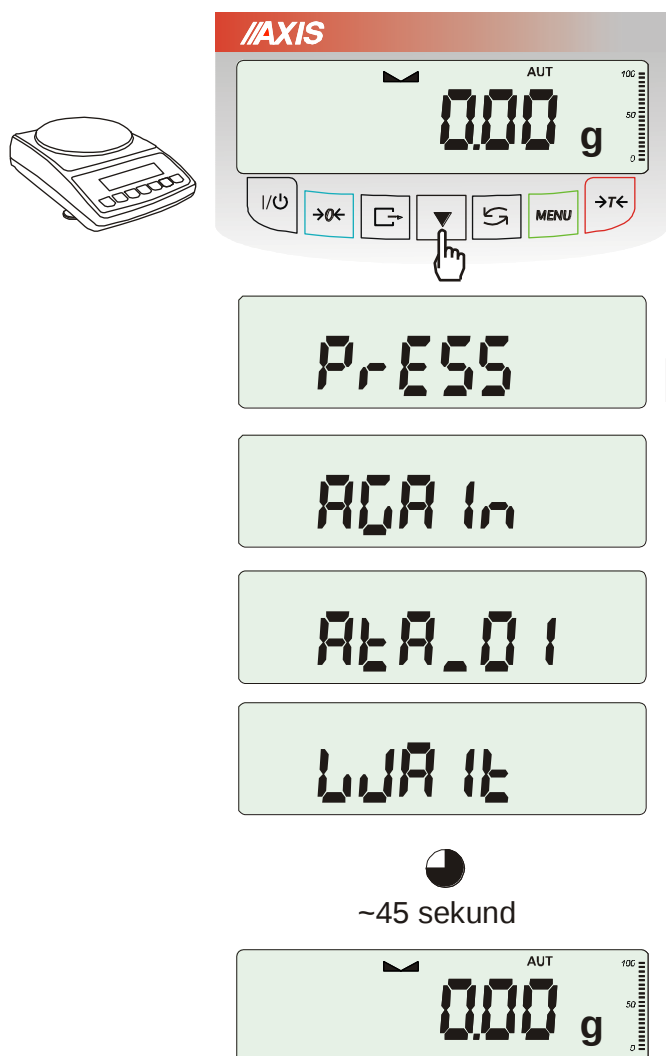

Empty the scale pan

Press the key ▼ twice (double-pressing the key is to avoid accidentally activating calibration).

During calibration, the internal weight is applied three times and the results are compared. Any discrepancies in the results are signaled with a message and weighing is suspended.

Do not perform any operations on the scale until the calibration process is complete. Any shocks or vibrations will disrupt the calibration process, extend its duration, and reduce the accuracy of the result.

If the user does not react for a longer period of time, automatic calibration will be initiated, preceded by a message *Aut_CAL*.

Correctly performed internal calibration ends with a zero reading when the scale pan is empty.

11. *Checking the balance*

In order to confirm correctness of the balance during its operation, before starting and after finishing every measurement series it is advised to check weighing accuracy. It can be done by weighing external calibration weight or other object with exactly known mass.

If exceeding of allowable measurement error is affirmed, the following things should be checked:

- if the balance stands stable and it is levelled,
- if the balance is exposed on rapid air blasts, vibrations, rapid temperature changes or air humidity,
- if the balance is not affected directly by heat source, electromagnetic radiation or magnetic field.

The cause of inaccuracy can be too low temperature of the balance as well, when it was unplugged from power supply. In this situation leave the balance switched on for several minutes in order to adjust its internal temperature.

If none of above causes of inaccuracy occurs, calibration with external weight should be performed to the balance. Recommended external calibration weight (to buy for additional charge) is given in technical data table. In order to calibrate the balance with external weight in legally verified balances verification seals should be removed and another legal verification should be performed. In this case it is recommended to contact authorized service centre.

Calibration with external weight is described in details in chapter 17.1.

12. Connection with a computer or a printer

The scale is equipped with RS232C, which can be used to connect external devices such as computer or a printer.

When cooperating with computer, the scale sends weighing result after initialize signal from computer or after pressing  key on the scale.

When cooperating with a printer data is send automatically after result stabilisation, but next transmission is possible after removing previously weighted sample.

When cooperating with label printer after pressing  key, the scale sends instructions set for the label printer. Label number 0001, hour, data (if the clock is installed and on) and nett weight. During transmission *LabEL* communicate is displayed.

The way of sending data and transmission parameters is set using *SErIAL* special function.

Set of send data is set using special function *PrInt*.

The following data can be send:

- Header (scale type, Max, d, e, serial number),
- Operator identification number,
- Successive printout number (measurement),
- Identification number or product bar code,
- Number of pcs (PCS function only),
- Single detail mass (PCS function only),
- Net weight,
- Tare (package mass),
- Gross weight,
- Total mass (Total function only).

If the scale is equipped with two serial joints *Print* function is set independently for both interfaces.


Computer must have a special program for cooperation with data from a scale.

Dedicated programs are also offered by AXIS.

Except RS232C joint, the scale can be equipped with USB or Wi-Fi interface.

Needed drivers and instructions are available on www.axis.pl.

Connecting cable *WK-1* (scale – computer / 9-pin interface):

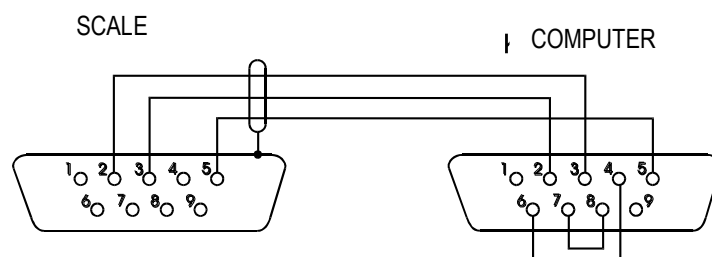

12.1 Detailed Long protocol description

In standard mode (Long protocol) transmission proceeds in the following way:

1. Communication parameters: 8 bits, 1 stop bit, no parity, baud rate 9600bps,
2. Available orders send from computer and balance answers:

- Readout of scale indication (corresponds to pressing  key

Computer→Scale: **S I** CR LF (53h 49h 0Dh 0Ah),

Scale→Computer: scale response according to description below (16 bytes):

Byte	1	-	sign „-“ or space
Byte	2	-	space
Byte	3÷4	-	digit or space
Byte	5÷9	-	digit, decimal point or space
Byte	10	-	digit
Byte	11	-	space
Byte	12	-	k, l, c, p or space
Byte	13	-	g, b, t, c or %
Byte	14	-	space
Byte	15	-	CR
Byte	16	-	LF

Attention:

Network number different than zero (*SERIAL / nr* function) changes scale working mode: communication with a computer is possible after logging the scale in with 02h scale number command. To log the scale out use 03h command.

For example: Using a program to test RS232 interface (program is available on www.axis.pl in computer programs section) for scale number 1 please write: \$0201 to log in, then SI, and write: \$03 to close communication.

- Asking for scale presence in system (testing scale connection with computer):
Computer→Scale: **S J** CR LF (53h 4Ah 0Dh 0Ah),
Scale→Computer: **M J** CR LF (4Dh 4Ah 0Dh 0Ah),
- Displaying a sign on scale display (text message from computer):
Computer→Scale: **S N** n n X X X X X X CR LF (53h 4Eh 0Dh 0Ah), nn-displaying time in seconds; XXXXXX- signs to display
Scale→Computer: **M N** CR LF (4Dh 4Eh 0Dh 0Ah),
- Scale tarring (calling $\rightarrow T \leftarrow$ key press) :
Computer→Scale: **S T** CR LF (53h 54h 0Dh 0Ah),
Scale→Computer: without response,
- Scale zeroing (calling $\rightarrow 0 \leftarrow$ key press):
Computer→ Scale: **S Z** CR LF (53h 5Ah 0Dh 0Ah),
Scale →Computer: without response,
- Scale turning on / off (calling $\uparrow \downarrow$ key press):
Computer→ Scale: **S S** CR LF (53h 53h 0Dh 0Ah),

Scale →Computer: without response,

- Entering to special function menu (calling *MENU* key press):
 Computer→ Scale: S F CR LF (53h 46h 0Dh 0Ah),
 Scale →Computer: without response,
- Setting low threshold value (option):
 Computer→ Scale: S L *D1...DN* CR LF (53h 4Ch *D1...DN* 0Dh 0Ah)
D1...DN – threshold value, maximum 8 characters („-“ – negative value, digits, dot – decimal separator), number of digits after dot should be the same as on scale display,
 Scale →Computer: without response,
- Example:
 - in order to set low threshold 1000g in scale B1.5 (d=0.5g) the following order should be sent:
 S L 1 0 0 0 . 0 CR LF (53h 4Ch 31h 30h 30h 30h 2Eh 30h 0Dh 0Ah),
 - in order to set low threshold 100kg in scale B150 (d=50g) the following order should be sent:
 S L 1 0 0 . 0 0 CR LF (53h 4Ch 31h 30h 30h 2Eh 30h 30h 0Dh 0Ah),,
- Setting high threshold value (option):
 Computer→ Scale: S H *D1...DN* CR LF (53h 48h *D1...DN* 0Dh 0Ah),
D1...DN – threshold value (see)
 Scale →Computer: without response.

12.2 EPL protocol description

Transmission parameters: 8 bits, 1 stop bit, no parity, baud rate 9600bps,


- After using  key in scale:
- Scale→Label printer : set of instruction in EPL-2 language that initialize label printing:

US	- Steering instruction
FR"0001"	- Label number define instruction
?	- Instruction that starts list of variable signs
mm:gg	- 5 signs: minutes:hour
rrrr.mm.dd	- 10 signs: year.month.day
masa	- 10 signs: scale indication+ mass unit
P1	- Steering instruction


Attention:

1. Except variable signs constant signs can also be inscribed e.g. factory name, product name and so on.
 2. In standard only one label pattern is possible to printout (number 0001). Using bigger amount of patterns (other label numbers) is possible thanks to *LABEL* special function.
 3. To achieve label printout, label printer must have inscribed label pattern (label pattern is created on computer and using computer it is saved to label printer memory). Label pattern is designed by ZEBRA DESIGNER program which is supplied together with label printer.
- Scales parameters and transmission protocol must correspond to label printer type.

13. Start-up

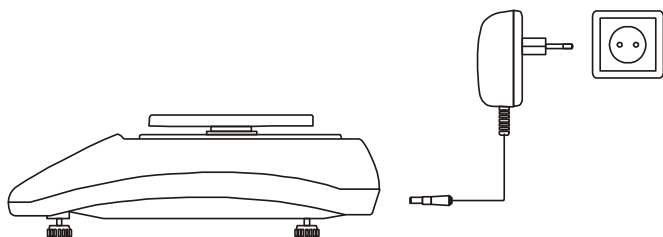


Plug feeder into ~230V power supply socket. When the pan is empty plug feeder output connector into 12V socket at back of the scale. Autotests and internal calibration will be performed.

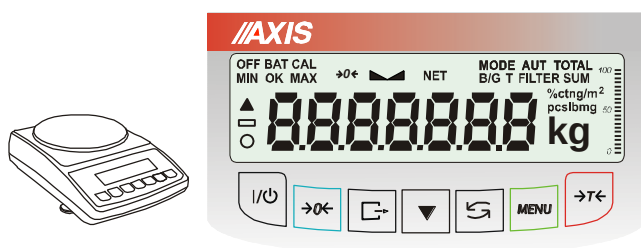

Autotest of balance display.



Showing scale program version.

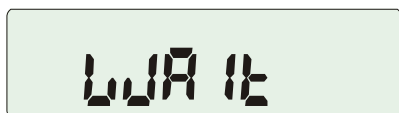

Internal calibration – 45 seconds (press ▼ key if You want to terminate calibration)




Ready to work.

Attention: *UnLOAD* communicate means that the balance is loaded or the transport securing elements were not clear away.



It is recommended that before you start measuring the internal temperature has stabilized weight. For this to happen, the weight should remain enabled for at least 2 hours. To maintain the accuracy of the weight it is not recommended to turn off the power.


14. Weighing with tare

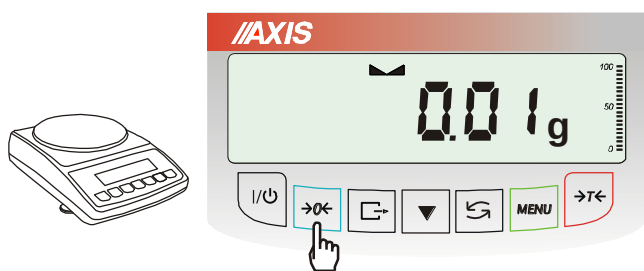

If the scale is not loaded and $\rightarrow 0 \leftarrow$ indicator doesn't indicate, press $\rightarrow 0 \leftarrow$ key.




Zero indication and $\rightarrow 0 \leftarrow$ indicator mean that the scale is ready to work.

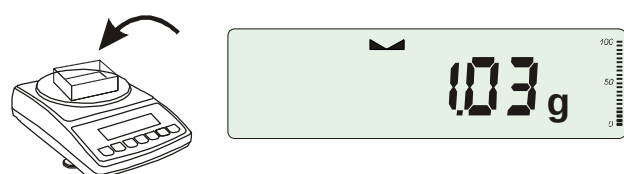


After putting container (package) tare the scale using $\rightarrow T \leftarrow$ key. NET indicator will show up.



Put on weighted object and readout net weight (NET indicator shows that scale indicates net weight).

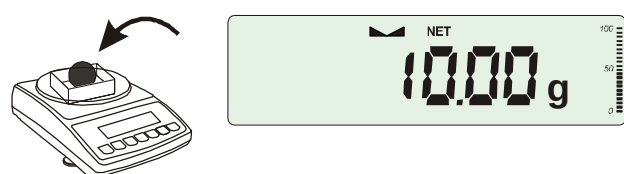


In order to readout gross weight press $\rightarrow G \leftarrow$ key (B/G indicator shows that scale indicates gross weight). Press again $\rightarrow G \leftarrow$ key in order to come back to net indications.

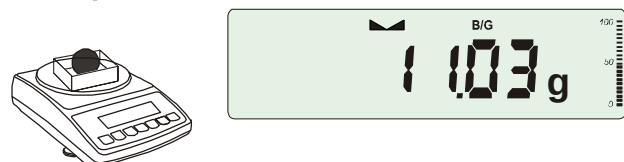

15. Scale menu

All scales except for basic metrological functions: weighing and taring, have many special functions and configuration options.

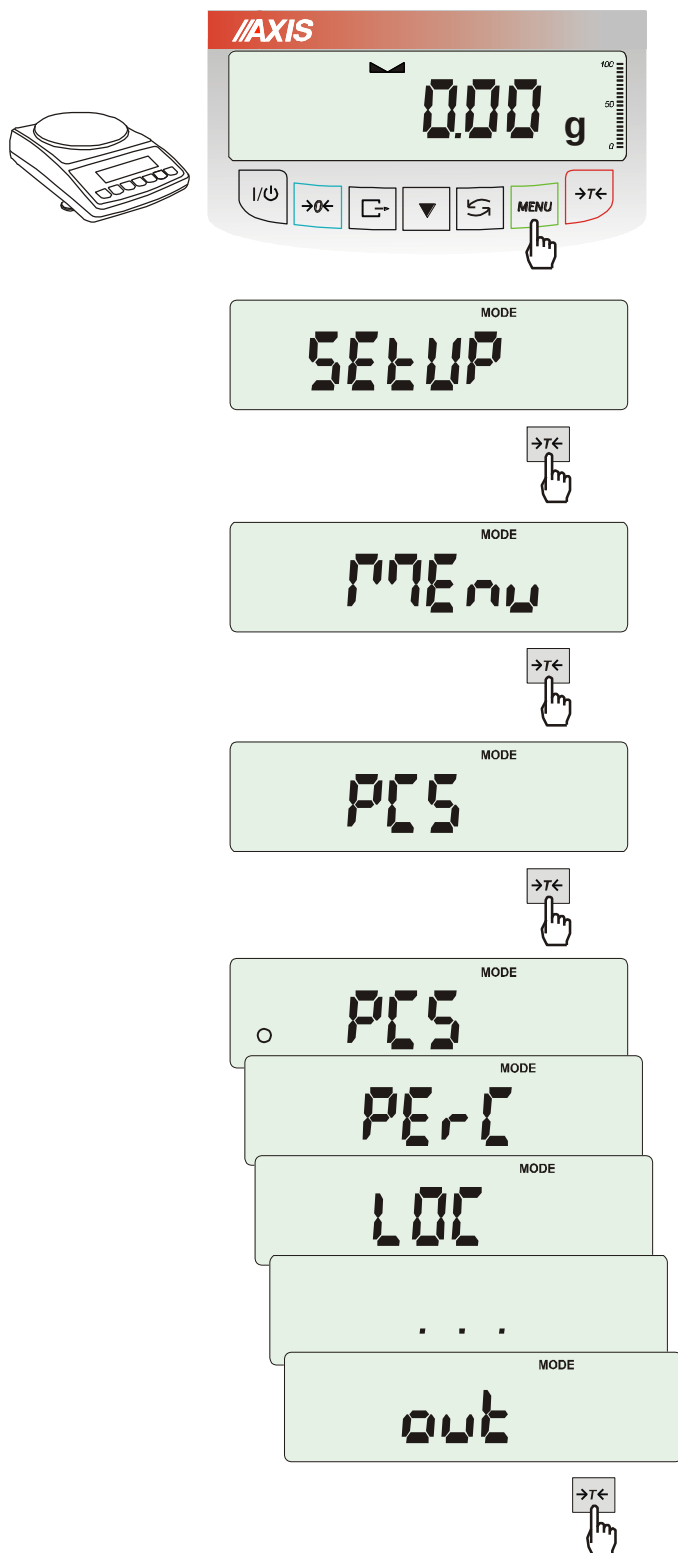

In order to ease using functions user can create his own (personalized) menu.

Creating personalized menu:

In „out of the box” scale after pressing *MENU* key only *SEtUP* option (it contains all configuration options) is available.




One of the configuration options is *Menu* that is used to create personalized menu.

To add a function to personalized menu press $\rightarrow T \leftarrow$ key when the function is indicating.

Chosen function is indicated with „o” sign on the left side of display.

After adding all necessary functions press *out* in order to come back to weighing mode. User now after pressing *MENu* key has access to selected earlier functions and to *SEtUP* option. *dEFAULT* option is used to set factory settings.

16. Menu navigation rules

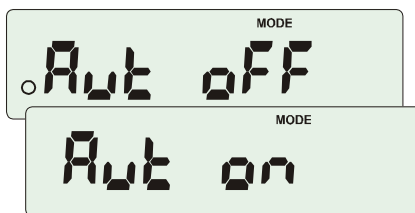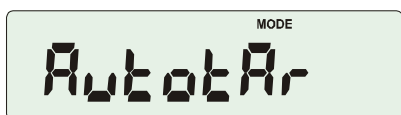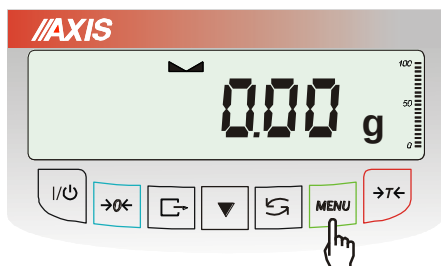

Choosing menu options:

Scale menu shows up after pressing *Menu* key. First menu position is displayed for about 10 seconds. After 10 seconds successive menu positions are displayed automatically.



Choosing menu position (option) is done by pressing $\rightarrow T \leftarrow$ key when it is displayed on the screen.

After choosing position (option) usually several options show up:

on – turning on selected option,

OFF - turning off,


out – out to menu.

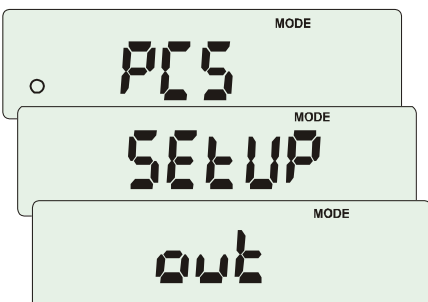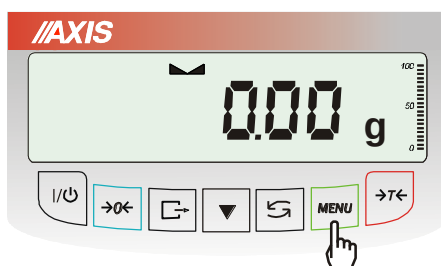


Accelerated working with menu:

First menu position is displayed for about 10s. User can change menu positions manually by pressing $\blacktriangledown$ key.



Immediate out to previous menu level is done by using *Menu* key.

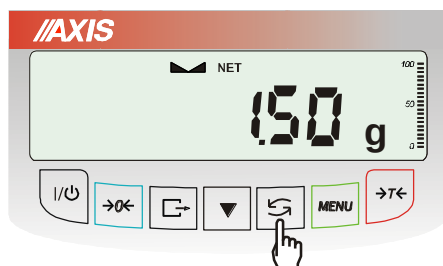


 key working method:

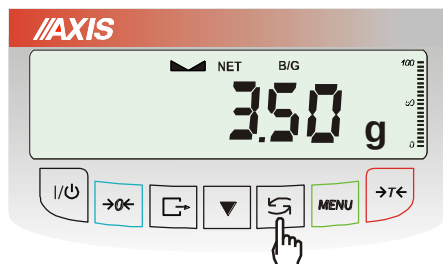


During standard weighing  key is used to switch between net and gross indication.

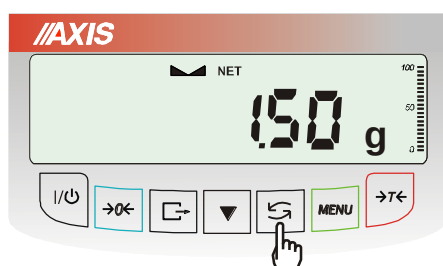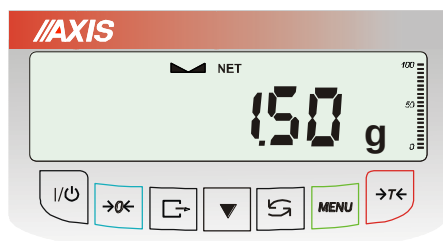


When special function e.g. PCS is turned on, using  key enables to go back to standard weighing mode.

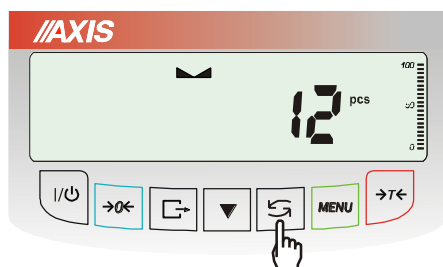

Sign „o” on the left side signalizes that special function is turned on and user can go back to function mode by pressing  key.

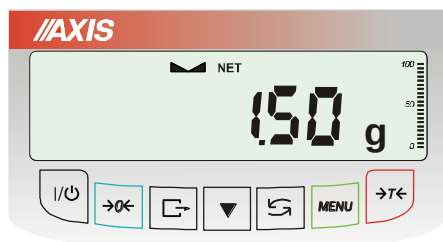

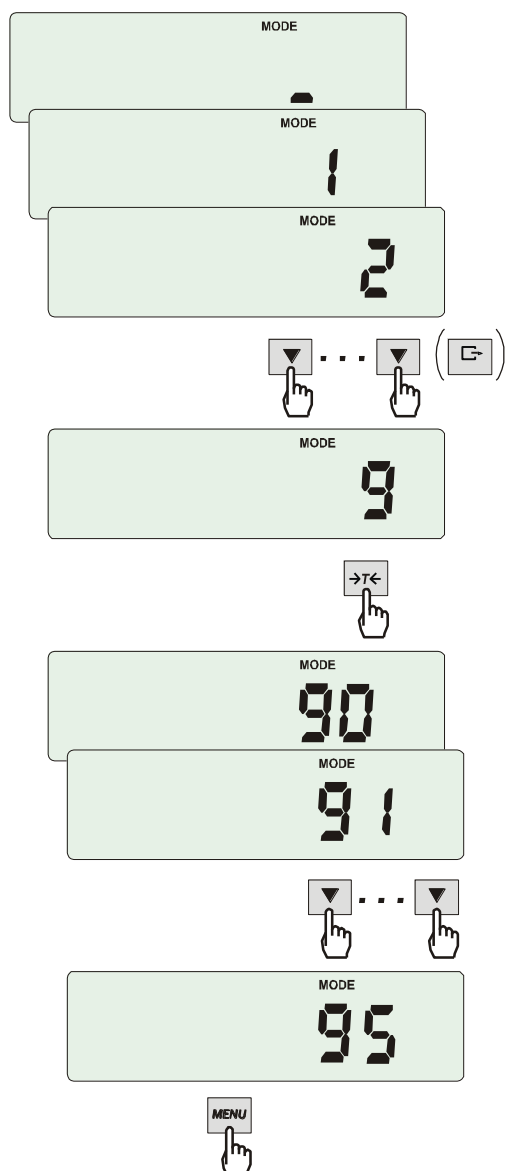

Inscribing numerical values:

Inscribing numerical values is needed in some special functions e.g. *tArE* function requires to inscribe tare values.


Keys:


 - decimal point,
 - next digit position,


 - increasing digit inscribed value (use also ) ,

MENU – end of inscribing.



Menu diagram:

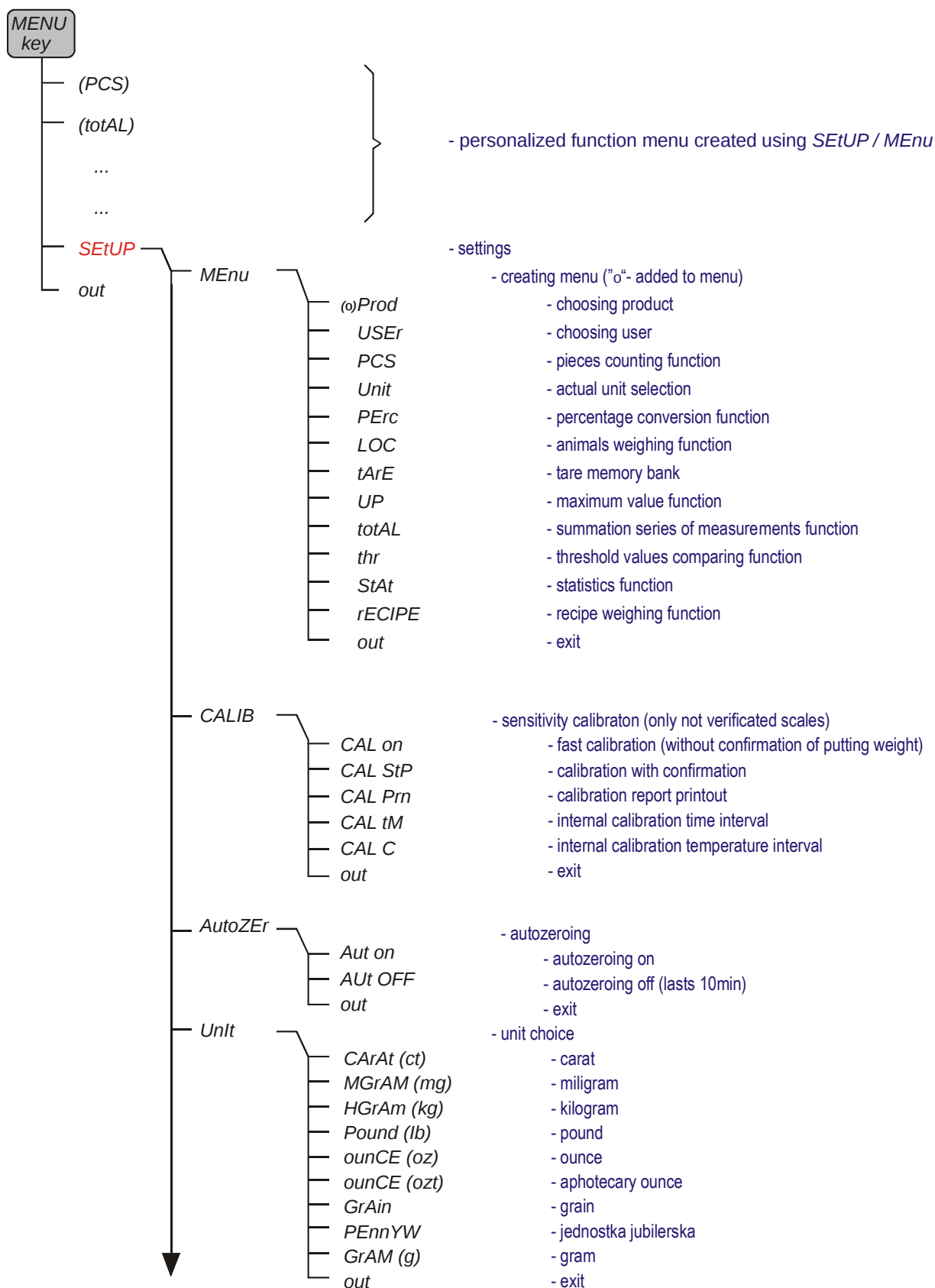

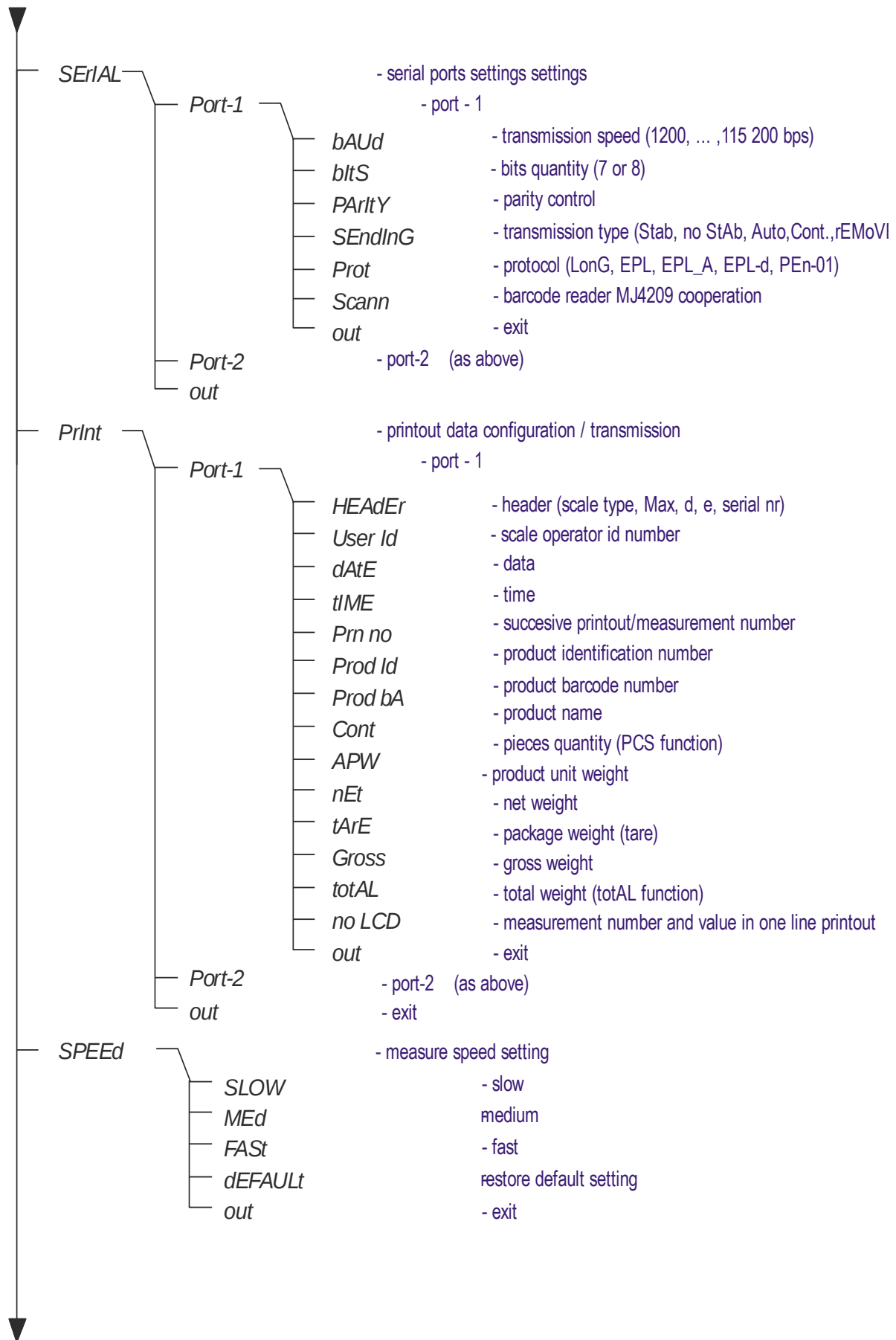

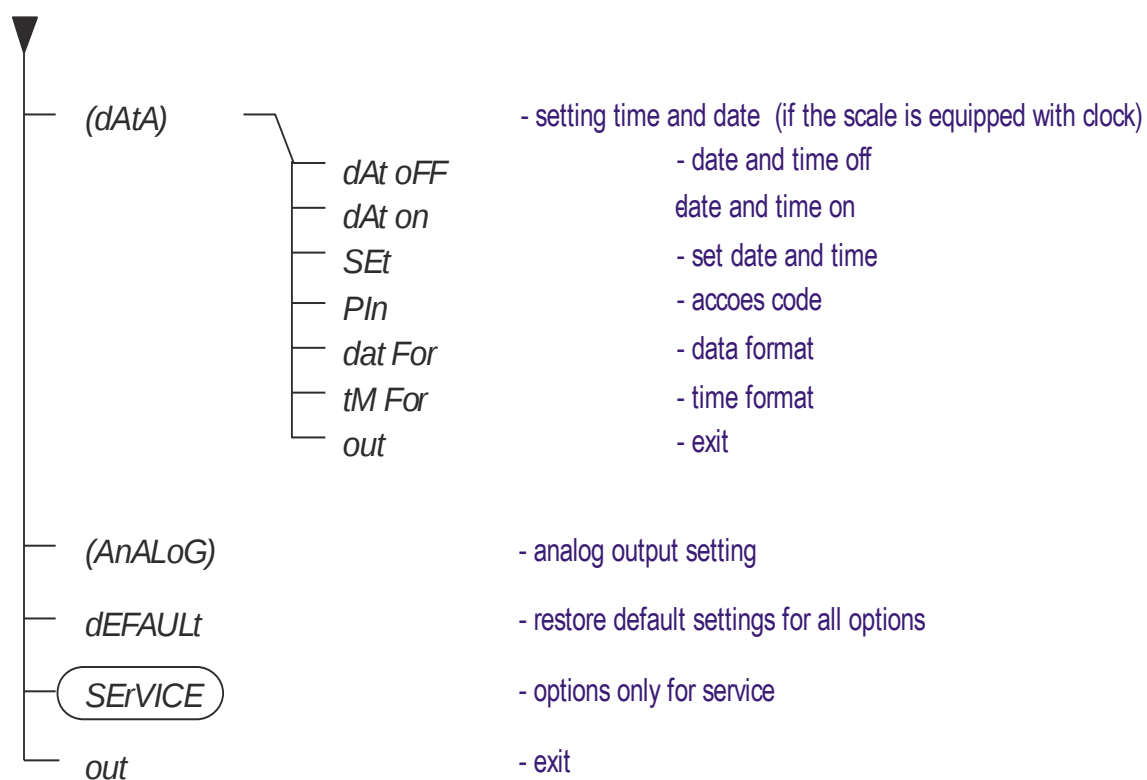

17. Scale setup (SEtUP)

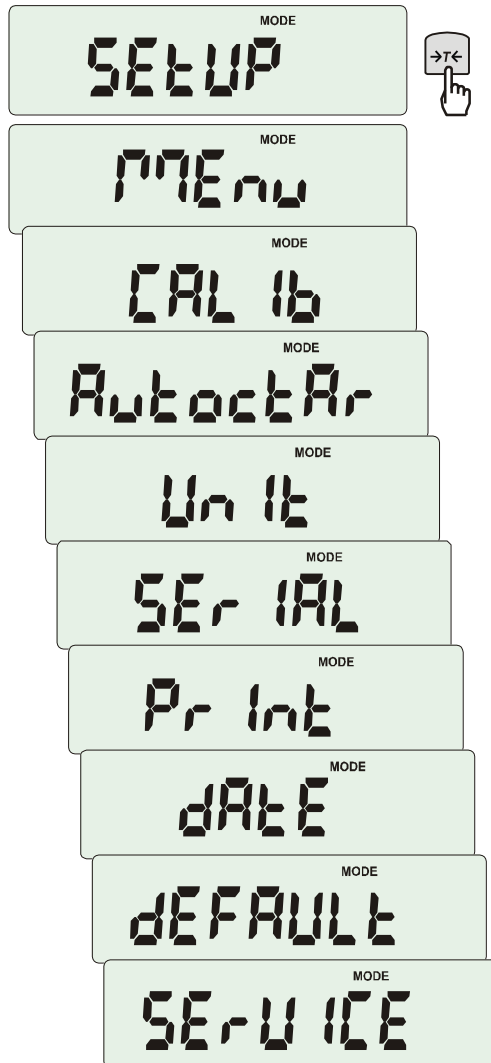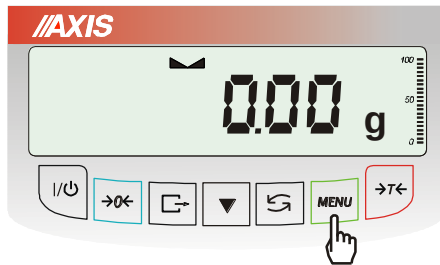

SEtUP contains all options used for setting scale work mode:


- ❑ MEnu – creating personalized user menu,
- ❑ CALIb – scale sensitivity calibration,
- ❑ AutoZer(ing) – self-maintaining zero indication (unloaded scale),
- ❑ UnIt – weight unit selection,
- ❑ SErIAL – setting serial ports,
- ❑ PrInt – transmission (printout) data selection,
- ❑ dAtE – inscribing actual date and time,
- ❑ dEFAULT – reset to factory settings,
- ❑ SErVICE – service menu (only for service).

17.1 External standard calibration and internal calibration options (CALib)

Calibration with external weight should be performed if balance accuracy after internal calibration is not satisfactory. Calibration weight stated in technical data table for the balance (or of better accuracy) with valid verification certificate should be used then.



Calibration of legally verified balance requires violating a mark used to protect an access to adjustment switch and results in losing legal verification. To renew legal verification of the balance, it is necessary to contact a service or notified body.

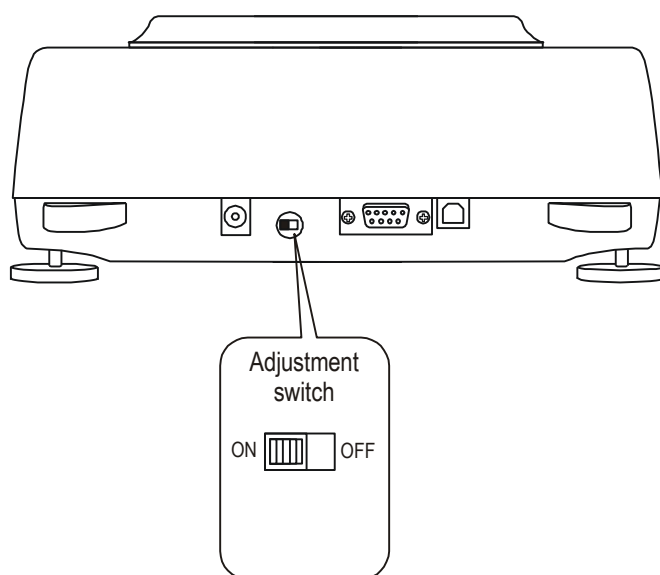


In balances comply with verification requirements performing calibration requires changing adjustment switch position, which is placed behind protecting mark (sticker) of a notified body. An access to the switch is possible only after removing the mark.

Before proceeding with calibration for balances comply with verification requirements, adjustment switch should be set to *ON* position using thin screwdriver (the balance will display the message *Pr ON*).

When calibration process, described on next page, is finished, the balance will display the message *Pr ON*. Adjustment switch should be set to *OFF* position using thin screwdriver (the balance will move to weighing).

Calibration with external weight:

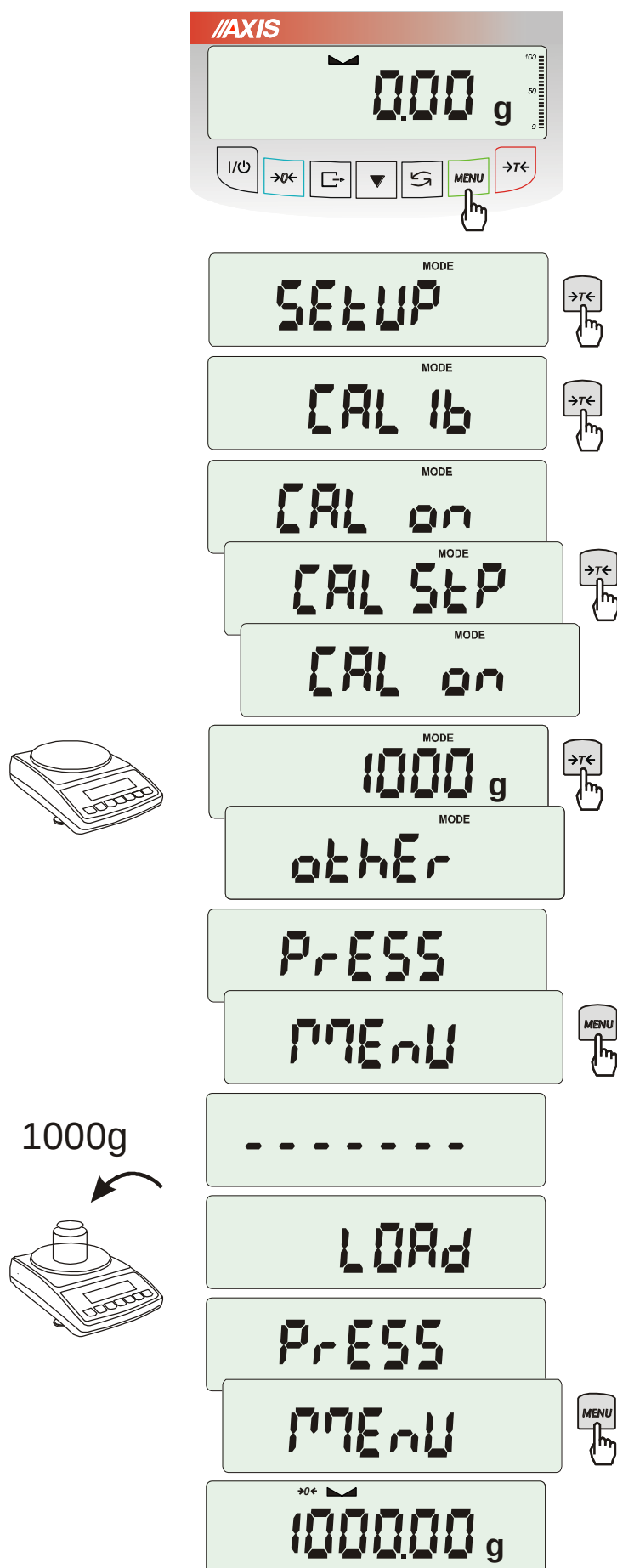

Press **MENU** key.

Press $\rightarrow T \leftarrow$ key when **CAL Ib** function appears.

The following options will be displayed:

- **CAL on** – calibration with external recommended standard of mass (see technical data).
- **CAL StP** – calibration with external weight, confirmation of successive steps - **MENU** key, **out** – leave without changes
- **CAL Prn** – calibration report,
- **CAL tM** – set time interval for internal calibration,
- **CAL °C** – set temperature difference for internal calibration,
- **out**.

Press $\rightarrow T \leftarrow$ key when **CAL StP** option appears (calibration in two steps).

Press $\rightarrow T \leftarrow$ key when weight value used for calibration is indicating or use **othEr** option and inscribe proper value (keys $\rightarrow 0 \leftarrow$, $\rightarrow \leftarrow$, $\rightarrow T \leftarrow$)

Press **MENU** and wait for writing zero to the scale (signaled by the message "-----").

When **LOAD** message appears put standard of mass on the pan. Press **MENU** key (**CAL on** doesn't need pressing **MENU** key).

Wait until internal calibration is finished and zero indication is displayed.

Note: Using **CAL on** instead of **CAL StP** frees you from pressing the **MENU** key twice

Internal calibration options:


Internal calibration options:

Internal calibration always takes place after turning on the scale (regardless of the options set), and can also be performed using a key ▼ or automatically:

- at specified time intervals (for verified scales – obligatory 2 hours),
- when the temperature changes by a given value (for verified scales – obligatory 2 °C).

In non-legally verified scales, the time and temperature values can be changed or turned off completely (the OFF option disables automatic internal calibration; only the key works ▼). When weighing is in progress, the start of internal calibration is delayed, but in legally verified scales, the delay time is limited to a few minutes.

Setting internal calibration options:

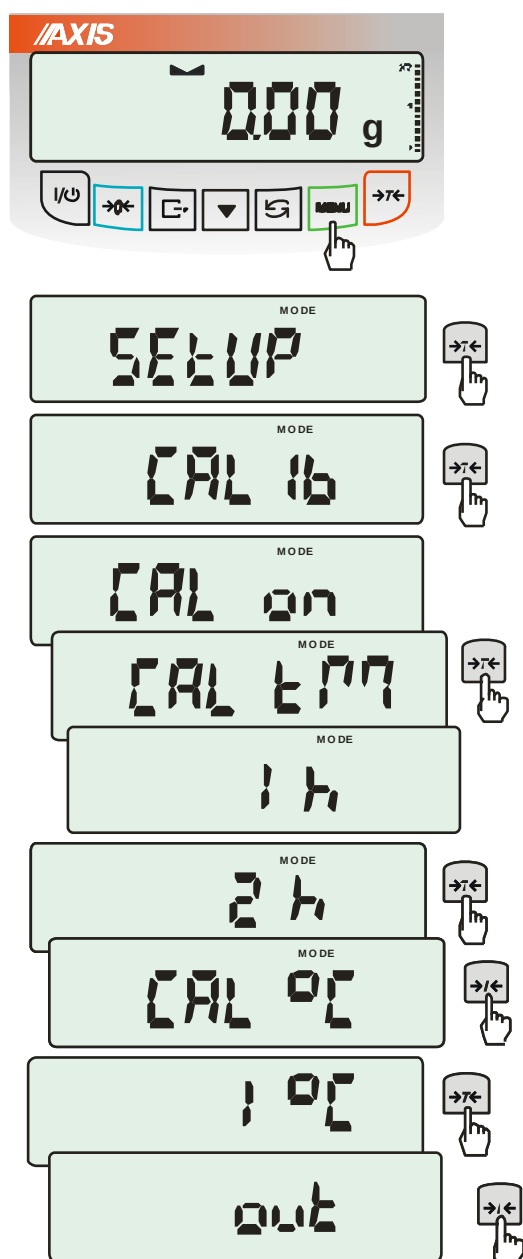

MENU key to call up the function menu and select the *CAL 1b* function by pressing the *T* ← key → when it is displayed.

The following options will appear:

- *CAL on* – performing calibration with an external weight standard
- *CAL Prn* – calibration report
- *CAL tM* – setting the time interval (OFF, 1h, 2h, 4h, 6h)
- *CAL °C* – setting the temperature difference (OFF, 1C °, 2C, 4C °)
- *out* – output without internal calibration

Press the →*T*← key while *CAL tM* is displayed. The internal calibration intervals will be displayed in sequence.

Select the desired value by pressing the →*T*← key.

Similarly, select *CAL °C* by pressing the key →*T*← and selecting the temperature difference value.

To exit, select *out*.

The form of ATA balance calibration report printout (option CAL Prn):

----- CALIBRATION REPORT -----

ATA2200 MAX=2200g e=0.1g d=0.01g
S/N : 1234
PROD.DATE: 2013-12-16
FIRM.VER.: ATA102 2015-01-12 AD7710 NTC

FACTORY EXT.LOAD : 2000.00 g
FACTORY INT.LOAD : 196.131 g
CALIBRATION NO. : 1
CALIBRATION DATE : 2015-01-22
CALIBRATION TEMP1: 30.346 'C
CURRENT EXT.LOAD : 2000.00 g
CURRENT INT.LOAD : 196.131 g
WEIGHT DIFFERENCE: 0.00 g


The form of ATZ balance calibration report printout (option CAL Prn):

----- CALIBRATION REPORT -----

ATZ2200 MAX=2200g e=0.1g d=0.01g
S/N : 1234
PROD.DATE: 2013-12-16
FIRM.VER.: ATA102 2015-01-12 AD7710 NTC

FACTORY EXT.LOAD : 2000.00 g
CALIBRATION NO. : 1
CALIBRATION DATE : 2015-01-22
CURRENT EXT.LOAD : 2000.00 g

17.2 Autozeroing function (AutotAr)

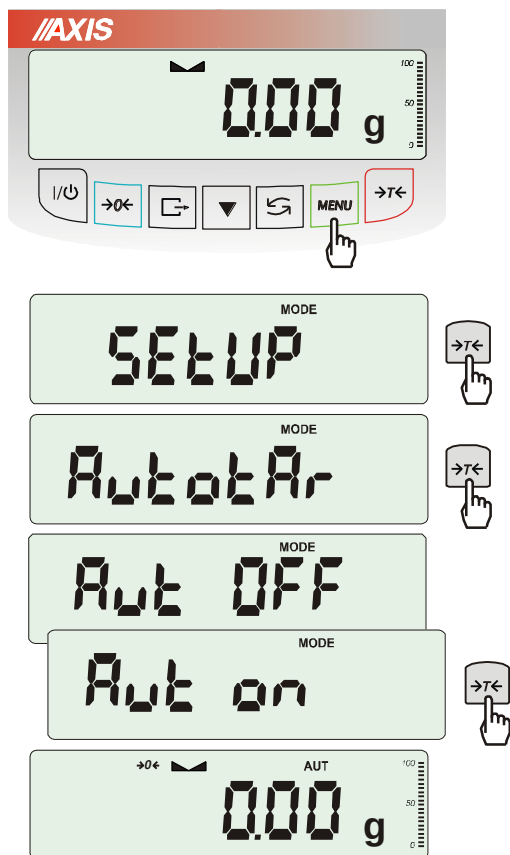

When the function is activated, the scale automatically ensures stable zero indication if the pan is empty or if zero indication was acquired by pressing $\rightarrow T \leftarrow$ key.



To turn on the function use **MENU** key and using $\rightarrow T \leftarrow$ key choose *AutotAr* and then *Aut on*

To leave the function press **MENU** key, then choose *AutotAr* and *Aut OFF*.

Note:

1. *AUT* sign occurs only in scales with LCD display.
2. In scales with active $\rightarrow 0 \leftarrow$ key function changes name into *AutoZE* (autozeroing) and also works when the zero reading is obtained by pressing the $\rightarrow T \leftarrow$ key.

17.3 Weight unit selection (Unit)

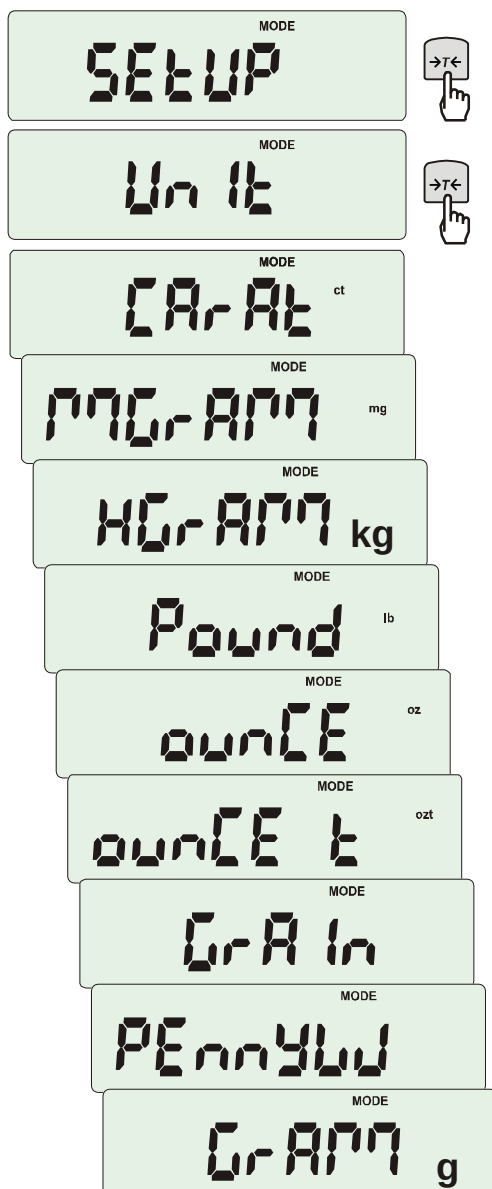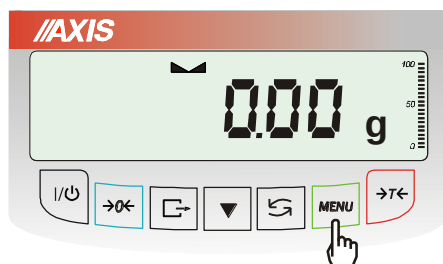

The function (in *SEtUP*) allows selecting weighing unit as default unit:

- CarAt (1 ct= 0,2 g) - carat,
- KgrAM (1kg=1000g) kilogram,
- Pound (1 lb=453,592374g) English pound,
- OunCE (1oz=28,349523g) - ounce,
- OunCEt (1ozt=31,1034763g) pharmaceutical ounce,
- GrAln (1gr=0,06479891g) - grain
- PennYW (1dwt=1,55517384g) jewellery mass unit,
- GrAM (1g) - gram.


The way of choosing carats as weighing unit is shown on the example.

Reading accuracy for individual units:

Unit	Readout	
	ATA220÷ATA520 ATZ220÷ATZ520	ATA1200÷ATA2200 ATZ1200÷ATZ2200
g	0,001 g	0,01 g
kg	0,000001 kg	0,00001 kg
mg	1 mg	10 mg
ct	0,005 ct	0,05 ct
lb	0,000 005 lb	0,000 05 lb
oz	0,000 05 oz	0,000 5 oz
ozt	0,000 05 ozt	0,000 5 ozt
gr	0,02 gr	0,2 gr
dwt	0,001 dwt	0,01 dwt

Select by pressing →T← when unit displayed.

17.4 Serial port parameters setting (SERIAL)

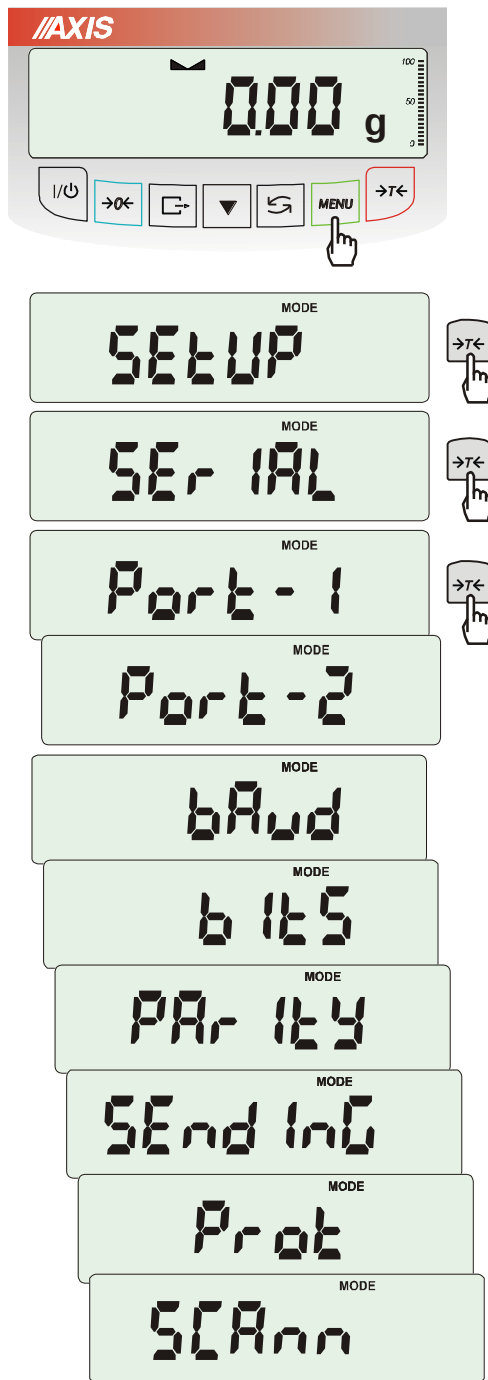


The function allows setting independently communication parameters of both of serial ports *Port-1* and *Port-2* (executed in RS232C, RS485, USB or LAN standard):

- transfer protocol (*Prot*):
 - Long* – cooperation with printer or computer,
 - EPL* – cooperation with label printer in normal mode (activates *LABEL* function),
 - EPL_A* – cooperation with label printer in automatic mode (activates *LABEL* function),
 - EPL_d* – cooperation with special label printers,
 - Pen-01* – cooperation with PEN-01,
- baud rate (*bAud*): (4800, 9600,115 200bps),
- number of bits in single char. (*bItS*): 7, 8,
- parity control (*PARItY*):
 - nonE* – no control
 - Odd* – nonparity
 - Even* – parity control,
- scale number in network (*nr*):
 - (if the scale doesn't work in network the number must be 0),
- transmission through serial interface (*SEndInG*):
 - StAb* – transmission after  key is used and result is stable,
 - noStAb* – transmission after  key is pressed without need of stabilisation,
 - Auto* – automatic transmission after load is put on and result is stable (*Auto*),
 - Cont* – continuous transmission, about 10 results per second (*Cont.*),
 - Remove* – transmission after putting off the weight.
- Default parameter values:
 - Long*, 9600 bps, 8 bits, *none*, *StAb*,
- *SCAnn* – cooperation with HD42A barcode readers.

In order to set needed parameters choose *SERIAL* function, select appropriate parameter and press *→T←* key when required option or parameter value is displayed. After setting the appropriate parameters and options, exit the function by selecting out.

In scales with an additional serial port appear *Port-1* and *Port-2*, for the independent setting of both ports.

17.5 Printout configuration (Print)

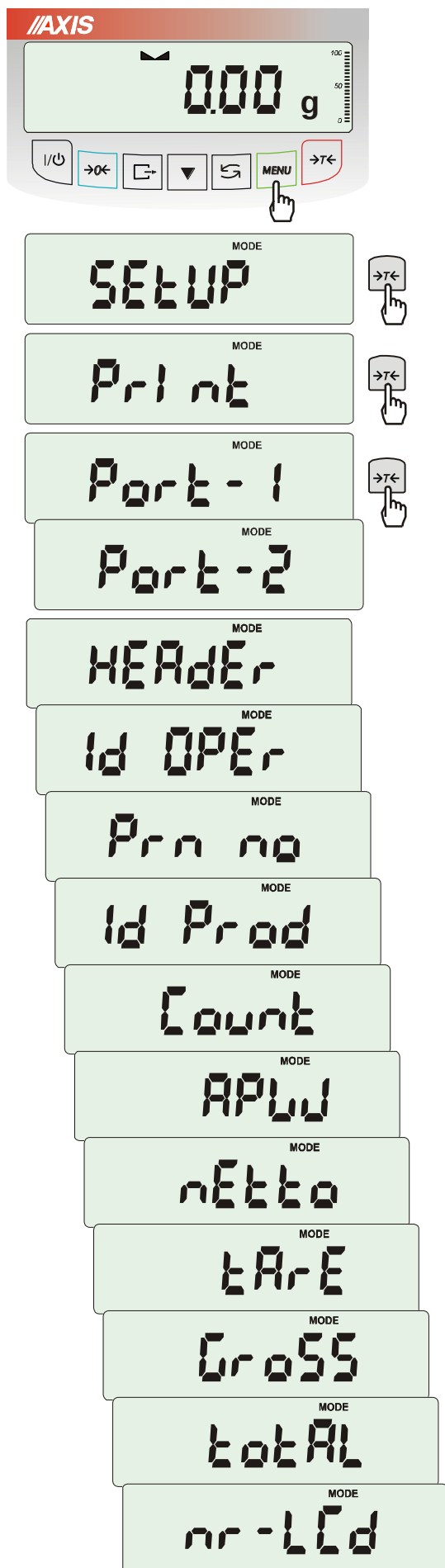

Function is used for printing additional information stored in scale memory, weighed product identification data and scale operator id. That information is inscribed using scale keys or scanner.

The function allows to switch on/off following positions on the printout:

- HEAdEr – header: name, model and scale number,
- USEr Id – scale user identification number,
- USEr nA – user name,
- Prn no – successive printout number (choose this option to zero counter),
- Prod Id – product number,
- Prod bA – product barcode (inscribed or scanned),
- Prod nA – product name,
- Count – counting result (PCS function),
- APW – unitary mass (PCS function),
- nEt – net mass
- tArE – current tare value,
- GroSS – gross mass,
- totAL – total mass (totAL function)

Attention:

If *Prod Id* or *USEr Id* is chosen, it is possible to inscribe quickly their new values (with omission of main menu).

In order to do that hold (about 3 seconds) *MENU* key and release it when *Prod Id* or *USEr Id* indicates. Inscribe new value using keys:

- 0← - increasing digit,
- T← - decimal point,
- T← - next digit,
- MENU* - end.

While inscribing *Prod id* user can use barcode reader connected to RS232C interface.

If the scale is equipped with two serial joints *Print* function is set independently for both interfaces.

Sample printout during normal weighing (all printout positions deactivated):

20.07 kg 20.04 kg 20.04 kg

Sample printout during normal weighing with clock option (all printout positions deactivated):


20.07 kg 2012-11-08 10:01 20.04 kg 2012-11-08 10:01 20.04 kg 2012-11-08 10:01

Sample printout during normal weighing (some printout positions activated):

ATA220	
MAX: 30kg e=d=0.01kg	
S/N :	
ID OPER.	: 000001
DATE	: 2012-11-08
TIME	: 12:26
NO	: 3
ID PROD.	: 01
COUNT	: 0 PCS
APW	: 0.000 g
NET	: 3.08 kg
TARE	: 0.00 kg
GROSS	: 3.08 kg
TOTAL	: 0.00 kg

17.6 Weighing speed selection (*SPEED*)

Option enables to change weighing speed, that enables better performance thanks to adaptation to environment conditions. In particular, if user chooses *FAST* speed setting then small weight portion batching for eg. in pharmacy is possible.

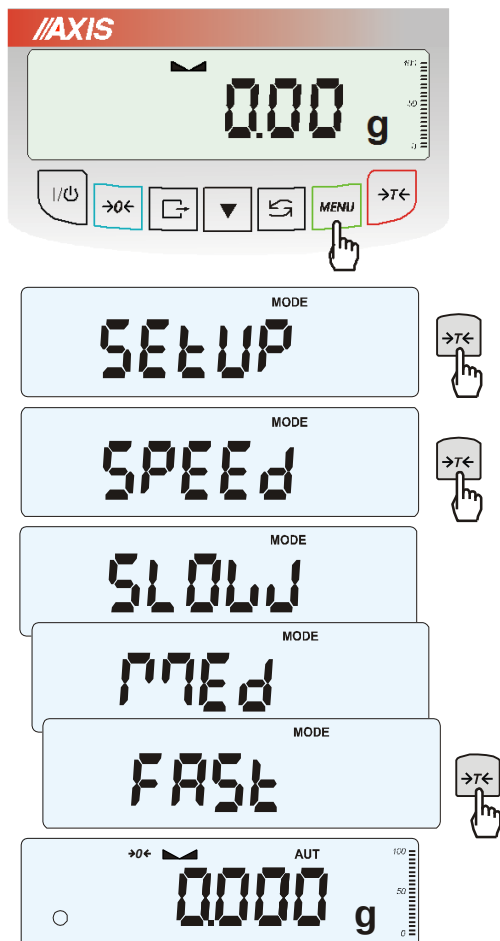

To turn on the function use *MENU* key and choose *SPEED* option by using **→T←** key, then select one of these options:

- *SLOW* – slow measurement (interferences),
- *MEd* – medium,
- *FAST* – fast (batching),
- *DEFAULT* – back to factory setting.

Attention:

When setting fast speed check if weighing results are stable. Otherwise use slower option.

18. Special functions description

All scales besides basic metrological functions: weighing and taring, have a set of special functions. Depending on meter type functions set differs. Below a list of functions available in standard ME-01 type meters:

- ❑ Add id number to product (*Prod*),
- ❑ Add id number to user (*USEr*),
- ❑ pieces counting function (*PCS*),
- ❑ change of mass unit (*UnIt*),
- ❑ percentage weighing function (*PErC*),
- ❑ selecting label number function (*LABEL*),
- ❑ weighing large animals function (*LOC*),
- ❑ entering tare function (*tArE*),
- ❑ maximum value indication function (*UP*)
- ❑ summing series of measurements (*totAL*)
- ❑ statistical calculations (*StAt*)
- ❑ recipe making (*rECIPE*)
- ❑ paperweight calculation function (*PAPEr*) – on demand

and functions that require additional equipment to be completely functional:

- options with the clock:
 - setting current date and time function (*dAtE*)
 - total weight function (*totAL*)
- options with the transoptors connectors (*WY᳚*):
 - checkweighing function (*thr*)

LabEL function is available in scales with *EPL* or *EPL-A* transmission protocol activates (go to *SetuP/SERIAL*).


18.1 Product and user identification (*Prod* and *USEr*)

The balance enables to inscribe product barcode and user identification number:


- *Prod bA* – product barcode,
- *USEr Id* – user identification number.

Balance product barcode and user readout together with external devices (e.g. printer, label printer and computer) enables to build simple identification and archivisation systems.


Inscribing multi-digit data without using e.g. computer keyboard is inconvenient and using barcode reader is beneficial.

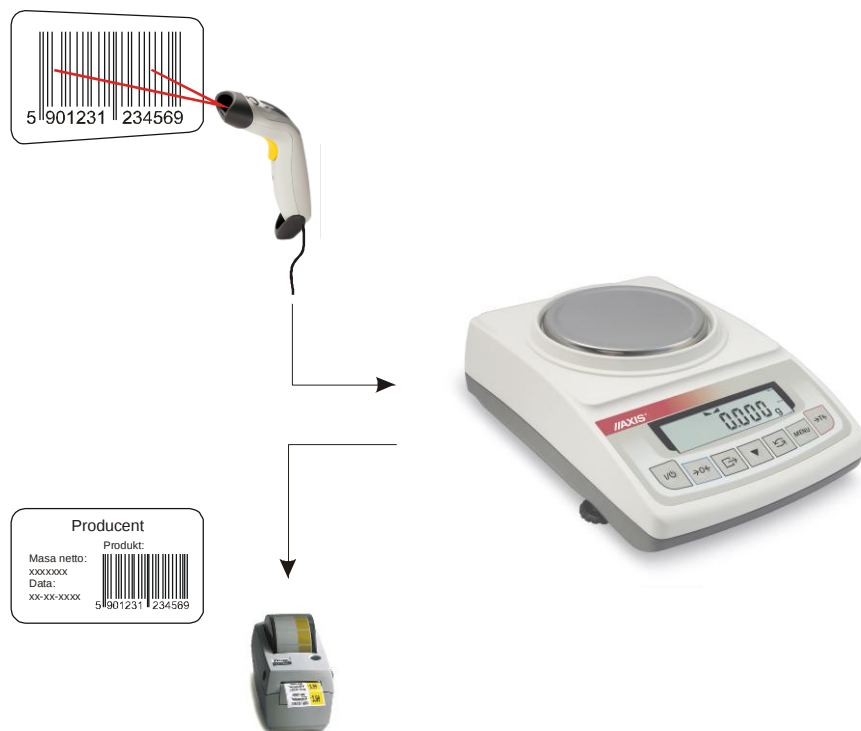

Prod and *USER* options enables inscribing single product and user data.

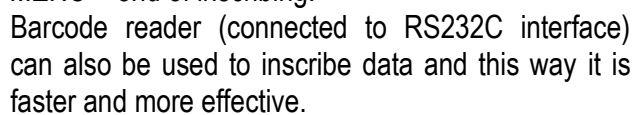

18.2 Pieces counting function (PCS)

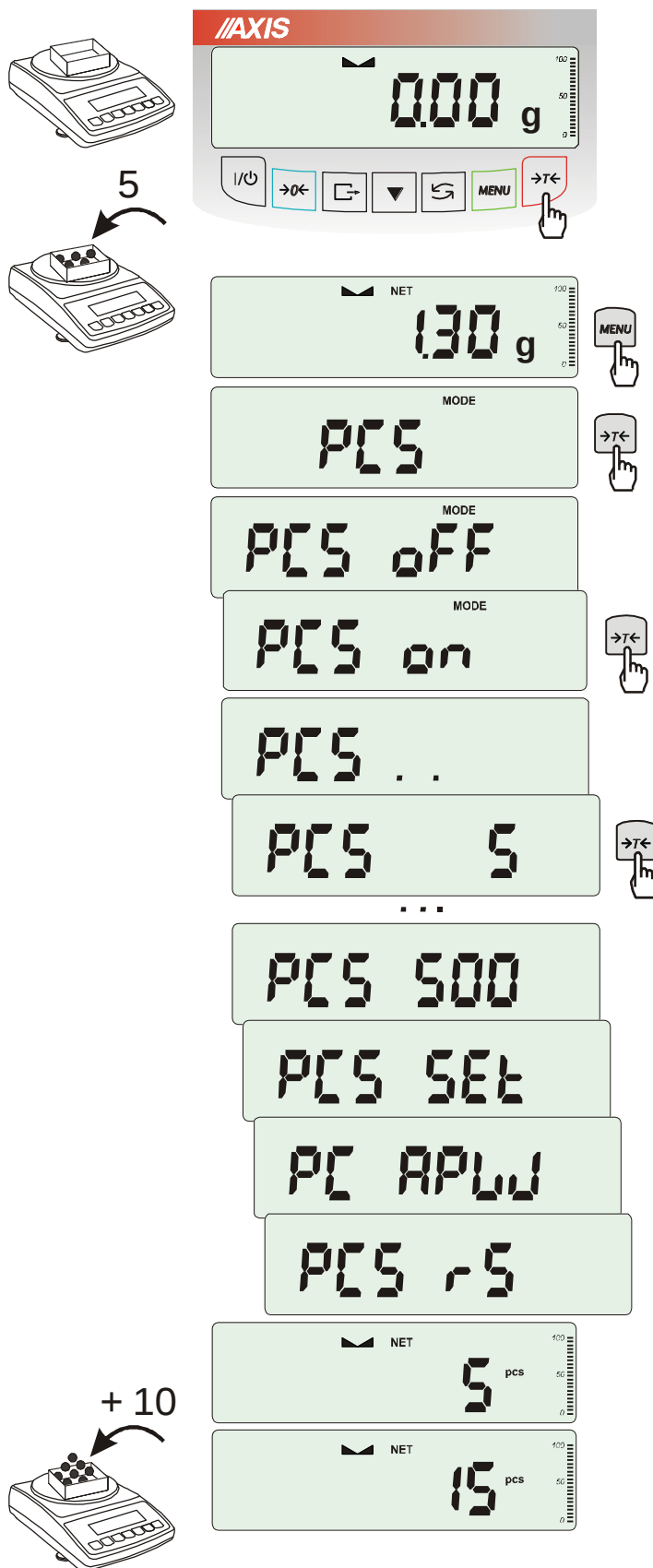

This function enables to count identical pieces, e.g. turnbuckles or buttons.

A measurement is performed in two phases:

- first phase - single piece weight calculation on the basis of defined pieces amount (5, 10, 20, 50, 100, 200 or 500 pieces),
- second phase – pieces counting.

First phase options:

- PCS . . – recalling of a value inserted earlier (this quantity must be inscribed earlier),
- PCS SEt – set any amount of pieces in a sample,
- PCS APW – set unitary mass directly,
- PCS rS – inserting number of details in a sample and receiving of their mass from other scale connected by RS-232C.


It is advised that single piece weight is not less than one reading unit and sample weight used in first phase is bigger than 100 reading units.

To leave function press MENU key and then using →T← key chose PCS and PCS OFF.

Note:

1. APW too LOW communicate signalises that a sample was not put on the pan or if single piece weight is less than one-tenth readout plot (counting is not possible).
2. APW LOW communicate signalizes that single piece weight is more than one-tenth but less than one readout plot. (counting possible but with bigger errors, result blinks).
3. In scales equipped with LED display pcs sign is replaced with "■".

18.3 Weight unit selection (UnIt)

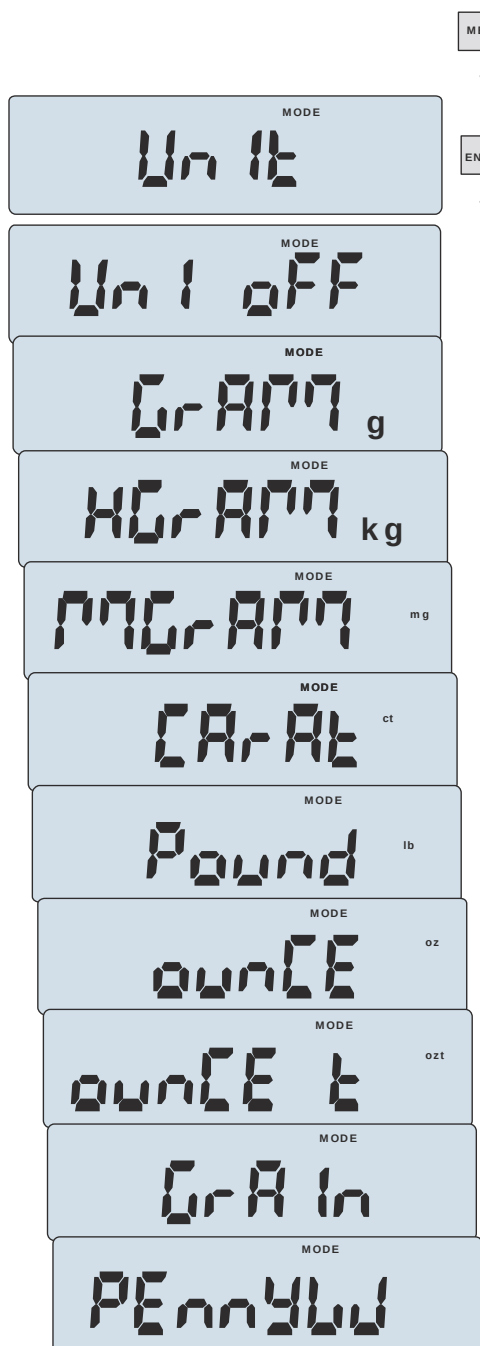

The function (in *Menu*) allows selecting current weighing (unit the power is turned off):


- GrAM (1g) – gram,
- kGrAM (1kg=1000g – kilogram),
- MGrAM (1mg=0,001g – miligram),
- CarAt (1 ct= 0,2 g) - carat,
- KgrAM (1kg=1000g) kilogram,
- Pound (1 lb=453,592374g) English pound,
- OunCE (1oz=28,349523g) - ounce,
- OunCEt (1ozt=31,1034763g) pharmaceutical ounce,
- GrAln (1gr=0,06479891g) - grain
- PennYW (1dwt=1,55517384g) jewellery mass unit,
- UnI oFF - return to default unit (set in *SEtUP*)

Reading accuracy for individual units:

Unit	Readout	
	ATA220÷ATA520 ATZ220÷ATZ520	ATA1200÷ATA2200 ATZ1200÷ATZ2200
g	0,001 g	0,01 g
kg	0,000001 kg	0,00001 kg
mg	1 mg	10 mg
ct	0,005 ct	0,05 ct
lb	0,000 005 lb	0,000 05 lb
oz	0,000 05 oz	0,000 5 oz
ozt	0,000 05 ozt	0,000 5 ozt
gr	0,02 gr	0,2 gr
dwt	0,001 dwt	0,01 dwt

Select by pressing →T← when unit displayed.

18.4 Percentage weighing function (PErC)

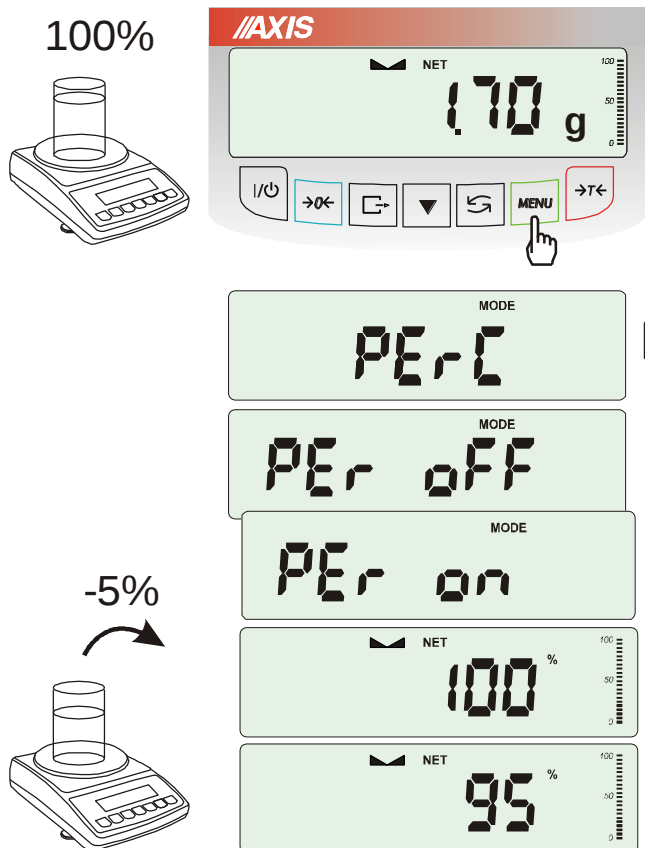

This function allows displaying weighing result in percents.

A measurement is performed in two phases:

- first phase – weighing a reference sample (100%),
- second phase – measuring specific sample as a percentage of the reference sample.

Weighing result is displayed in different format, depending on the weight value of reference sample.

The function has the following options:


- *PEr OFF* – disable the function,
- *PEr on* – set current scale indication as 100% and activate percentage weighing,
- *out* – exit without changing settings.

Note:

PEr Err message informs that reference 100% mass is less than 0,5*Min or was not defined.

18.5 Label choosing function (LAbEL)

This function is used in scale with *ELTRON* (*SErIAL* function) data protocol. This protocol enables label printout with actual scale indication and chosen data from *PrInt* special function (variable data), for example date and time. Other data, for example company address, product name, barcode can appear on label as a constant text. Label patterns with number (4 digit) used by user should be saved in scale memory according to printer manual. Label pattern choice is made by inscribing label number using *LAbEL* function.

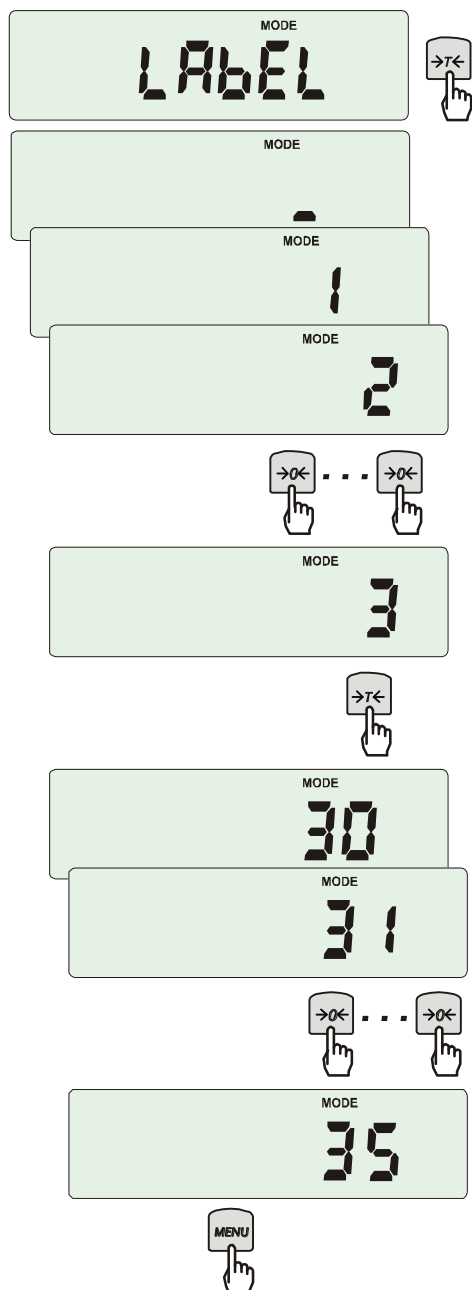

Press *MENU* button.

When *LAbEL* is displayed press →T← key.

Actual label number will show.

To enter new label number press →T← key, to exit function without number change press *MENU*.

To inscribe label number use keys:

▼ - digit increase,

→T← - next digit,

MENU – end.


After entering label number, putting load and pressing  key will cause sending data to label printer.

Data format sent to label printer (label nr 35, language EPL-2):

US	(55 53 0D 0A)
FR"0035"	(46 52 22 30 30 33 35 22 0D 0A)
?	(3F 0D 0A)
00:00	(30 30 3A 30 30 0D 0A)
2000.00.00	(32 30 30 30 2E 30 30 2E 30 30 0D 0A)
10 g	(20 20 20 20 20 31 30 20 20 67 0D 0A)
P1	(50 31 0D 0A)

18.6 Weighing animals function (LOC)

The function allows weighing animal moving on the scale.


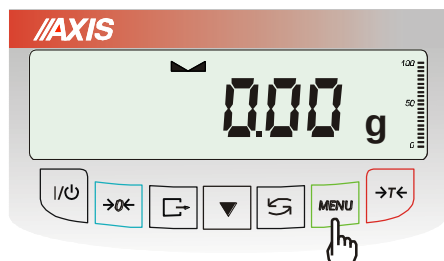


Press **MENU** key.

When **LOC** function is displayed press **→T←** key.

The following options appear on display successively:

- **LOC oFF** – leave the function,
- **LOC on** – automatic weighing after loading the scale,
- **LOC Pm** – the measurement initiated manually by pressing **↵** key.

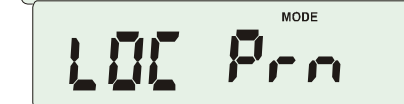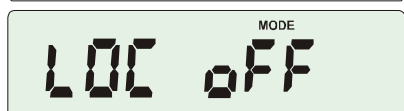


When **LOC on** is displayed press **→T←** key.

Tare the scale using **→T←** key if necessary and place the animal on the pan.

Wait until the weighing result is averaged – scale display blinks. Then scale will show stable (averaged) result and will send it through serial port.

The result remains on display for about 30 second.

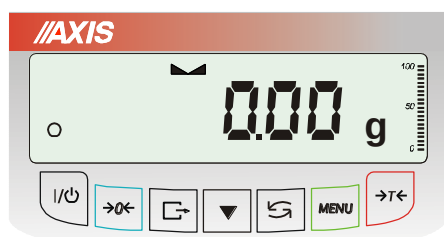


Important notes:

1. The loads lower than Min value are not averaged.
2. In case when putting animal on scale takes more than 5s it is suggested to choose **LOC Pm** option (measurement started manually by pressing **↵** key).

18.7 Constant tare memory function (tArE)

This function enables to measure gross weight of a sample placed in a container of a known weight value (stored in the memory) and to display calculated net weight of the sample. Tare value is recalled from the memory with $\rightarrow 0 \leftarrow$ or $\rightarrow T \leftarrow$ key when the pan is empty. Tare value may be entered using keypad or by putting container on the pan.

Inscribing tare value to memory:

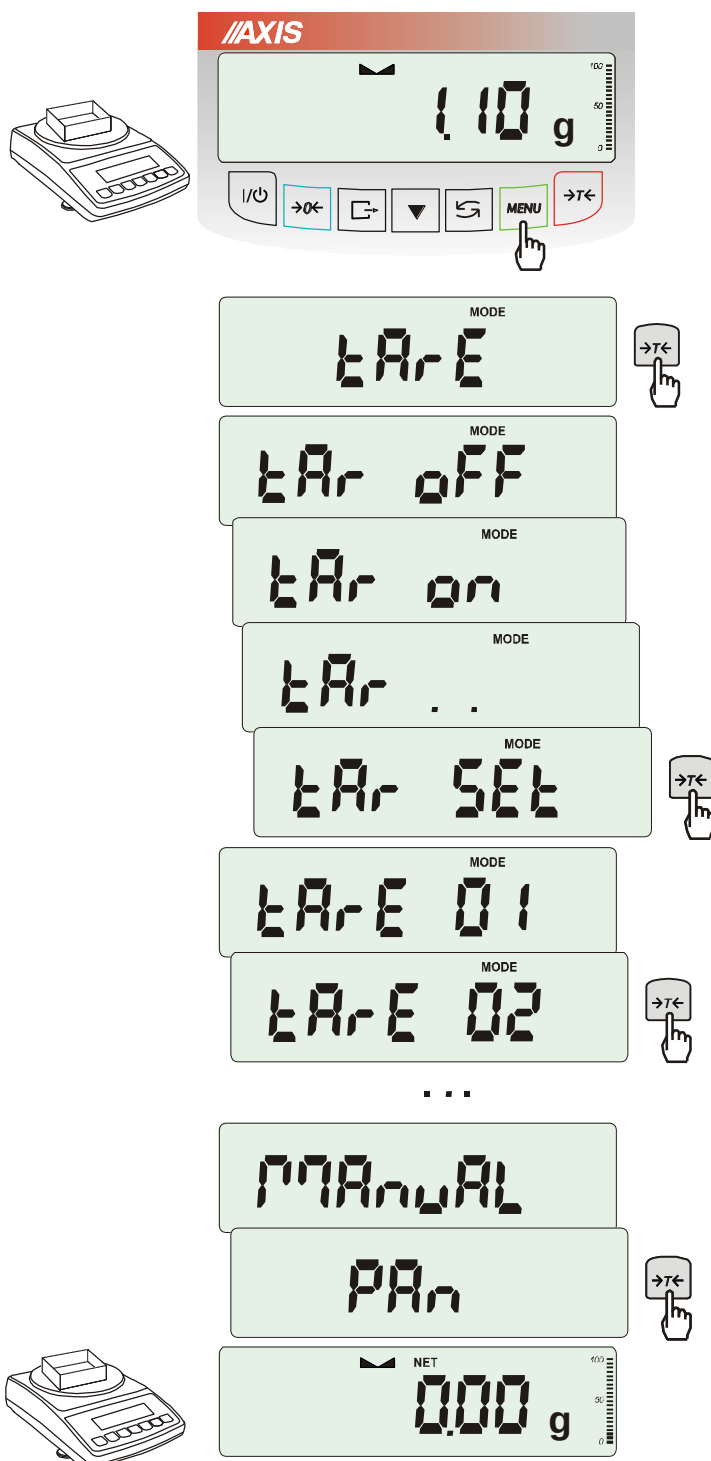

After pressing **MENU** key and choosing **tArE** function using $\rightarrow T \leftarrow$ key, the following options are available:

- **tAr OFF** – leave the function,
- **tAr on** – activate the function with the previous tare value,
- **tAr ..** – sample tare value from the pan,
- **tAr SET** – enter tare value with keys: $\rightarrow 0 \leftarrow$, $\rightarrow T \leftarrow$ and **MENU**
- **out** – printout a setting value of tare.

Press $\rightarrow T \leftarrow$ key when **tAr SET** is displayed.

By pressing $\rightarrow T \leftarrow$ key choose proper memory cell where tare will be stored: **tAr 01, 02, ... , 10**.

Choose inscribing method :


- **MANUAL** – inscribing using keys: ∇ , $\rightarrow T \leftarrow$ and **MENU**,
- **Pan** – inscribing mass value that is on the pan.

After storing tare, the scale starts working with inscribed tare value.

Note:

Tare value is stored in memory also after unplugging the scale from the mains.

Weighing with constant tare:

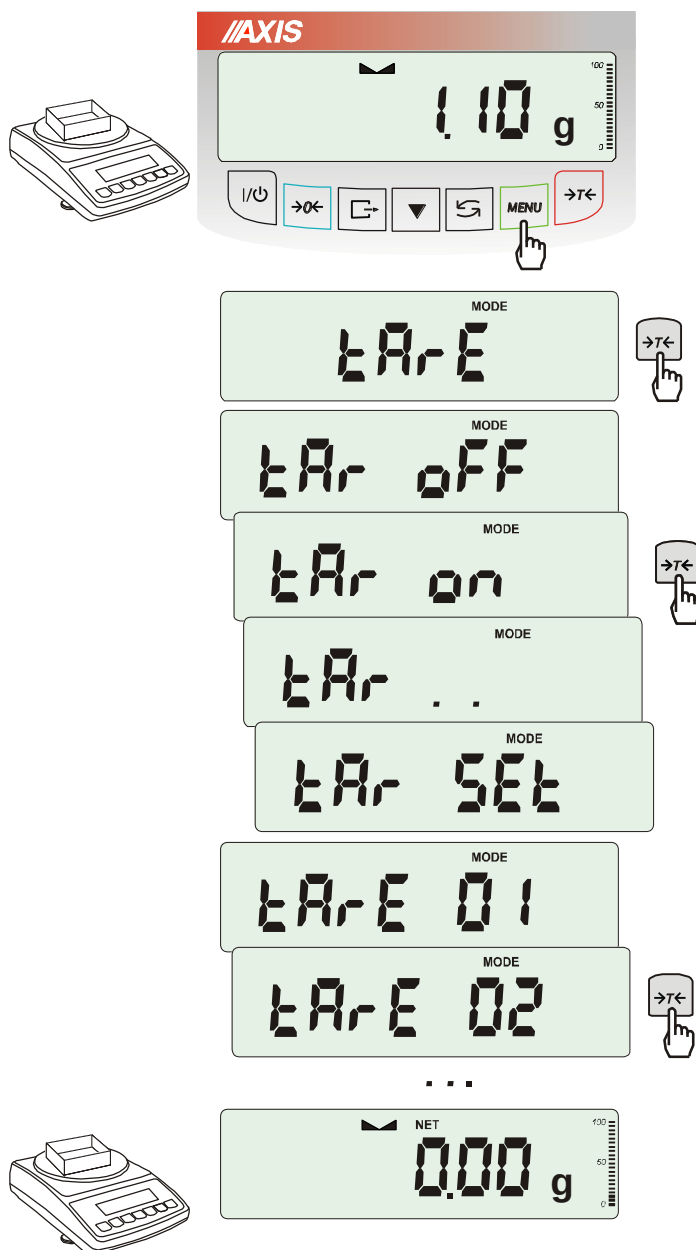

In order to use tare value that is located in memory, choose from menu *tArE* function and then *tAr on* option.

A list of memory cells will show up:

tAr 01, 02, ... , 10.


Cells with inscribed value are marked with "o" sign on the left side, active value marked with "▲".

Choose proper memory cell using **→T←** key.

tArE function is activated with chosen tare value. Moreover the scale will indicate nett weight (weight on the pan minus tare values). Using **→T←** key (or **→0←**, while empty pan) causes scale zeroing and then subtraction of recalled tare. Minus indication will show up.

18.8 Maximum value indication function (UP)

This function allows holding maximum (or minimum) value that is indicating at the moment.

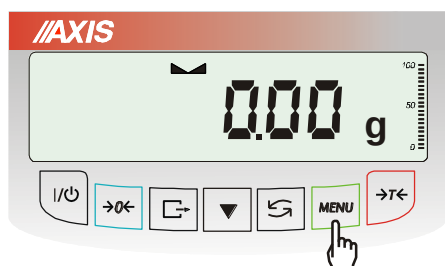

Before measurement scale should be tared.


Function has following options:

-UP oFF – function off,


-HIGH – holding maximum value,

-LOW – holding minimum value.

Pressing →T← key will cause result zeroing.

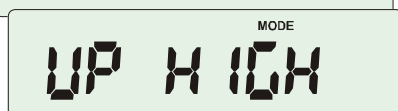









...

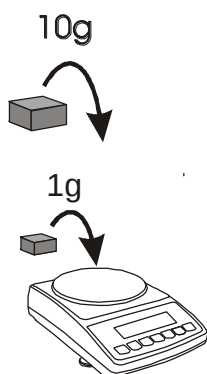

Note:

Autozeroing function and the stabilisation indicator are deactivated when UP function is running.

18.9 Total weight function (totAL)

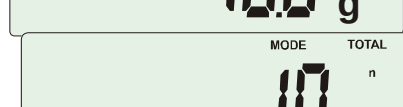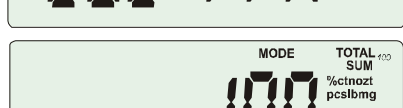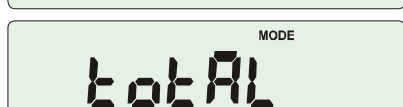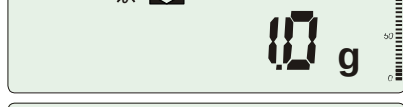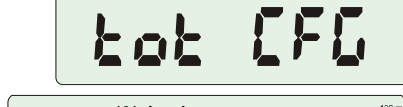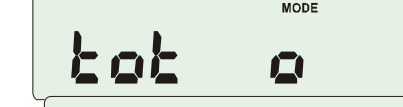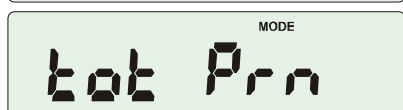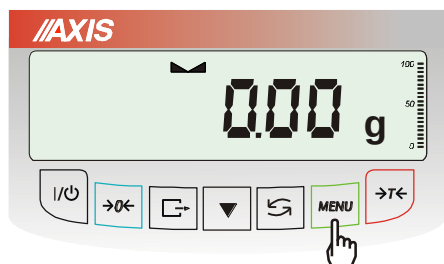

The function allows calculating total weight for series of measurements, which can be greater than scale capacity. It allows calculating total weight as well as average value.

Press **MENU** key.

When **totAL** is displayed press **→T←** key.

The following options will appear successively:

- **tot Prn** - report printout without clearing total register,
- **tot oFF** - clearing total register, report printout and leaving the function,
- **tot □** - working with receipt printout after each measurement,
- **tot -** working without receipt printout,
- **tot CFG** - saving measurement mode (using **↵** key: *Manual*, after taking off the load : *auto*).

Press **→T←** key when **tot □** is displayed.


Perform measurement series by pressing **↵** key for storing results into total register.

In order to print and display results enter the function by choosing **totAL** and **tot Prn** option from menu.

The results are displayed in the following sequence:

- total weight (SUM),
 - number of registered measurements (n),
 - average value (=),
- regarding that moving to display successive result is performed after pressing **↵** key.

In order to go back to total weighing without zeroing total register press **↵** key several times.

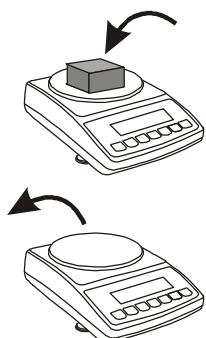

To leave the function with clearing total register, select *totAL* function from menu and choose *tot oFF* option. Scale prints the communicate informing about clearing registers.

The form of receipt after each measurement:

Date: ...	Time. ...
measurement no	weight
measurement no	weight

Report form:

Date: ...	Time. ...
TOTAL WEIGHT	=
NUMBER OF SAMPLES	=
AVERAGE VALUE	=

Note:

When the scale doesn't have an internal clock, Date and Time do not appear on printout.

Maximum number of measurements is 99 999.

Maximum total load 99 999 000d.

The weighing unit of the total value from the register (Total) is the same as the weighing unit stated on the keypad or is 1000 times greater, what is signalled by "o" indicator at the left of the display.

If the registered value is too big to be displayed, "E" communicate appears on the display. If the number of series is too high and cannot be displayed, "Err1" communicate appears on the display

18.10 Checkweighing function (thr)

This function allows comparing weighing result with two programmed reference values: lower and upper threshold. Comparison result is signalled with indicators (MIN, OK, MAX) and sound signal generated when threshold values are exceeded.

If comparison result is:

- smaller than zero threshold – no signal,
- smaller than lower threshold – the scale signals MIN (yellow colour),
- between threshold values - the scale signals OK (green colour, with the short sound signal),
- greater than upper threshold - the scale signals MAX (red colour, long sound signal).

The checkweighing results can be use to control:


- optical indicator (*Indication mode*),
- batching devices (*Batching mode*).

-











Standard scale is set for cooperation with optical indicator.

On outputs P1-P3 (*Relays socket*) short-circuit states appear as result of comparison scale indication with threshold values.

On the chart below output states are shown during increasing load on the scale for both working modes:

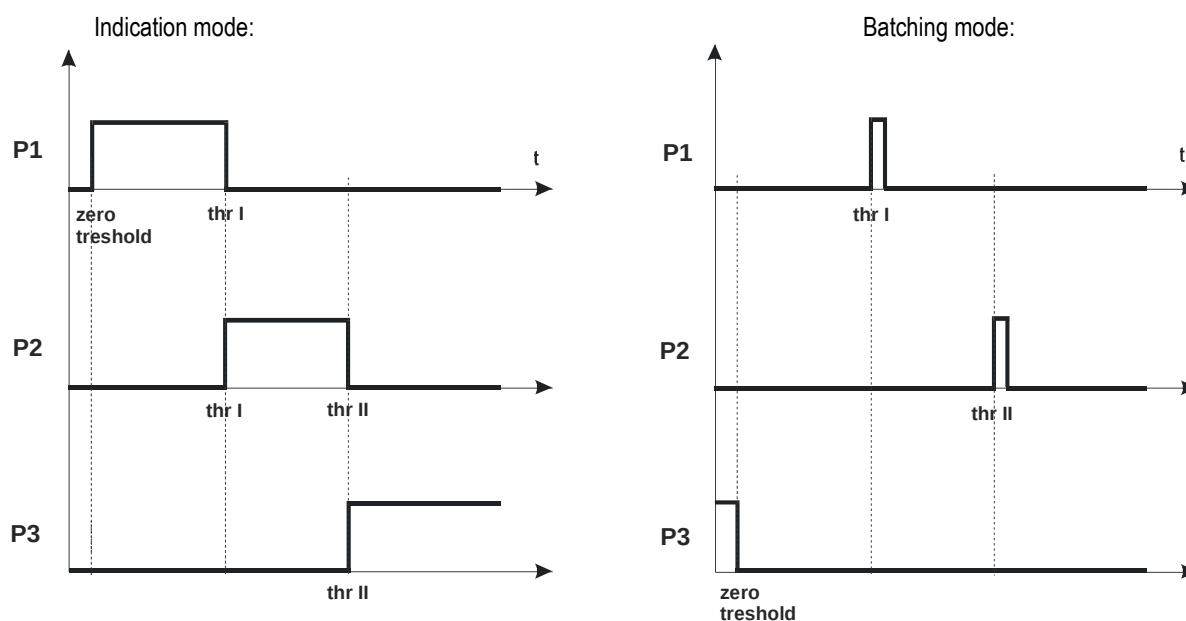

In *Batching* mode on P1 (thr I) and P2 (thr II) outputs short-circuit impulses appears for time of 0,5s. On P3 (zero) output short-circuit state appears when indication does not exceed threshold value signalling zero load.

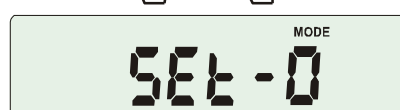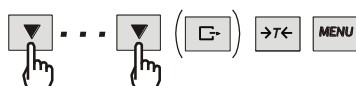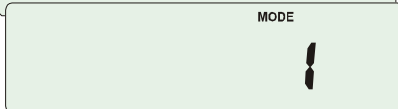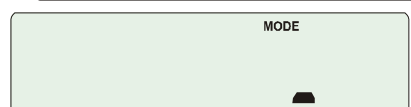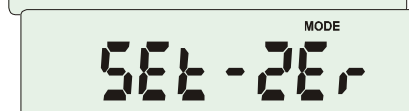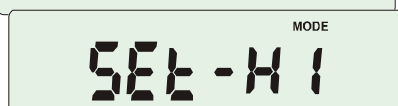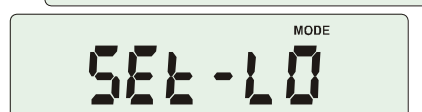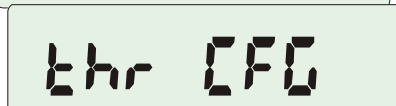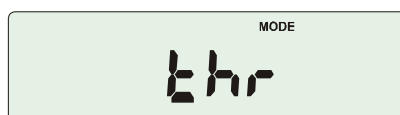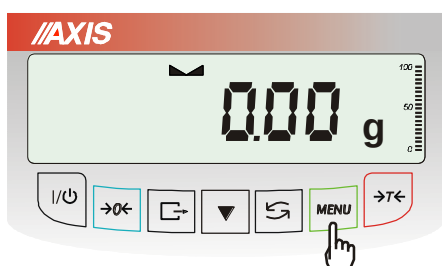Operation sequence:

Press **MENU** key and choose *thr* pressing **→T←** key.

The following options are displayed successively:

- *thr off* – deactivate the function,
- *thr on* – activate the function,
- *thr Prn* – check last threshold values (press **↵** key several times),
- *thr CFG* – choose *Relays* socket mode:
 IMPULS - *Batching* mode
 SIGNAL – *Indication* mode.

Choose *thr-on* option using **→T←** key. The following options for entering thresholds are displayed:

- *SEt-LO* - set lower threshold value,
- *SEt-HI* - set upper threshold value,
- *SEt-ZEr* - set zero signalisation threshold.

Using **→T←** key select *SEt-LO* option.

Set lower threshold value using the following keys:


- ▼** - digit increase,
- ↵** - decimal point,
- T←** - move to next digit,
- MENU** - finish.

Then select *SEt-HI* option and enter upper threshold value.


Choosing *out* option will cause starting work with signalisation of exceeding thresholds and zero.

To change *Relays* socket mode use *thr CFG* option. Default option is *Indication*.

To leave the function, press **MENU** key and then choose *thr* and *thr off* options.

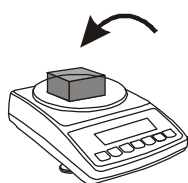


Relays socket (WY□):

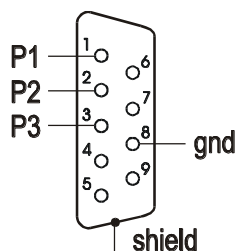

Relays output (WY□) is the open collector transistor output with load capacity 25mA / 24V.

Relays connection diagram:

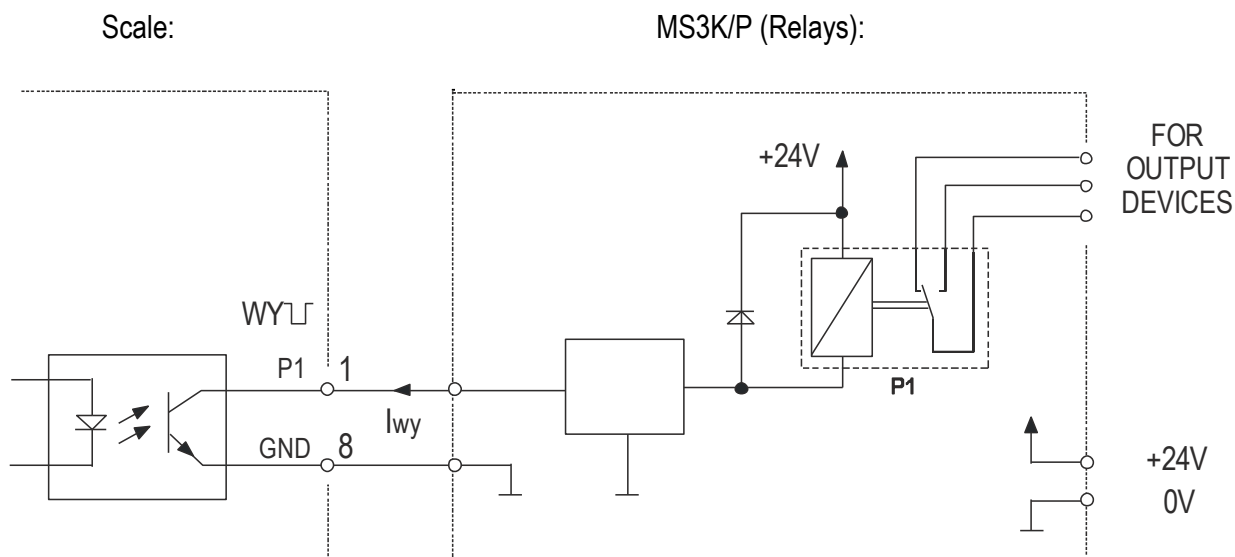


It is advised to use MS3K/P electronic board (sold separately), consisting of RM96P transmitters, with DC24V input voltage and AC250V, 3A output (for output devices).

Notes:

1. After switching the scale on, both thresholds are set to maximum values.
2. When setting upper threshold value, pay attention that its value is not below lower threshold value.
3. Setting lower and upper threshold value is possible after sending appropriate orders from computer, what is described in 12.1 Detailed Long protocol description.

18.11 Setting date and time function (dAtE)

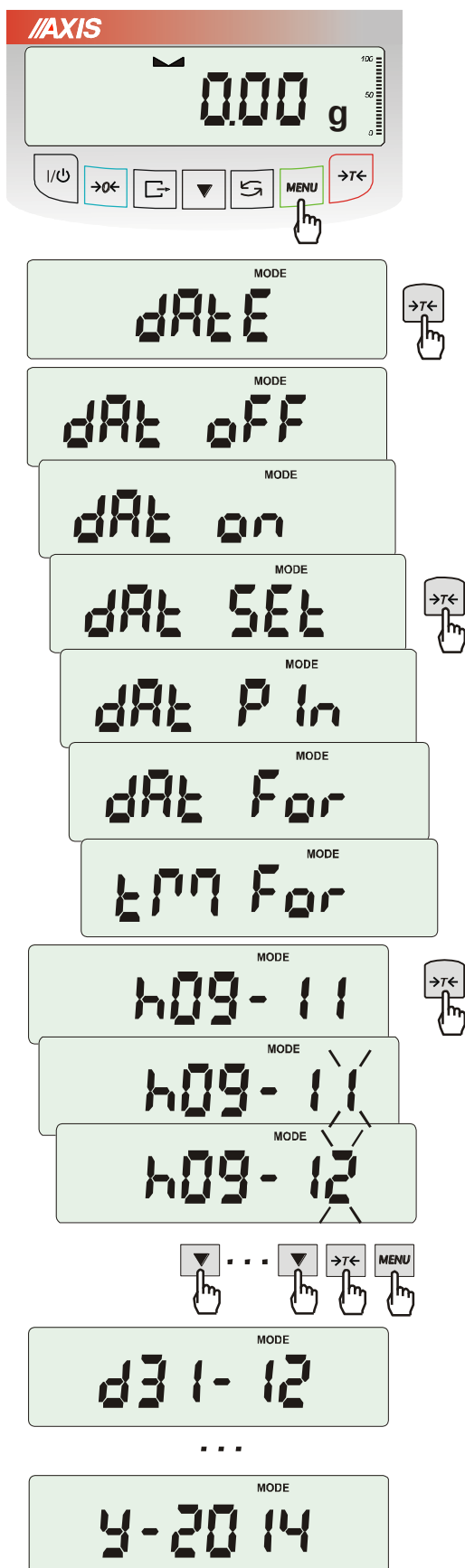

The function allows setting current date and time of scale internal clock and mode of its use.

The function has the following options:

- *dAt OFF* – deactivate date and time during printout of current weighing result,
- *dAt on* – activate date and time during printout of current indication ($\rightarrow$ key),
- *dAt SEt* - change current date and time,
- *dAt PIn* – data and time secure password (to prevent from changing date and time by unauthorized personnel),
- *dAt For* – data printout in *USA* or *EU* format
- *tM For* – time printout in 24h or 12h format.

The example at the left presents how to set current date and time using *dAt SEt* option.

After setting proper date and time activate it with *dAt on* option.

UE: rrrr-mm-dd gg:mm

USA: mm-dd-rrrr gg:mm AM/PM

(gg – hours, mm – minutes, AM – before noon, PM – after noon, mm – month, dd – day, rrrr – year).

Attention: Insccribing non-zero *PIN* value causes showing *PIN* sign during next date and time changing and inscribing 4 digit code is necessary. (using keys ∇ , $\rightarrow T \leftarrow$ and *MENU*).

18.12 Statistical calculations function (StAt)

This function evaluates from series of measurements (max 1000) statistical parameters of weighting process.

Adding successively measurements to register is automatic and it occur after the scale is loaded and its indications stabilize.





After each loading printout is made with: number of measurements, result, date and time (if clock is installed and the function is activated).

For the obtained measurements series the scale evaluates:

- n -number of samples
- sum x -sum of all samples $sum_x = \sum x_n$
- $\bar{x}$ -average value (sum x)/n
- min -minimal value from n samples
- max -maximal value from n samples
- max-min -maximal value minus minima value
- S -standard deviation $S = \sqrt{\frac{1}{(n-1)} \sum_n (x_n - \bar{x})^2}$
- srel -variance factor $srel = \frac{S}{\bar{x}}$

Statistical calculations results can be printed.

Order of operations:

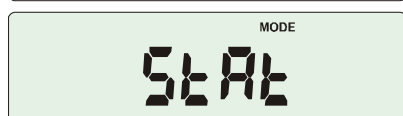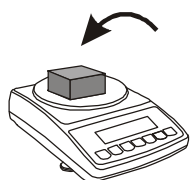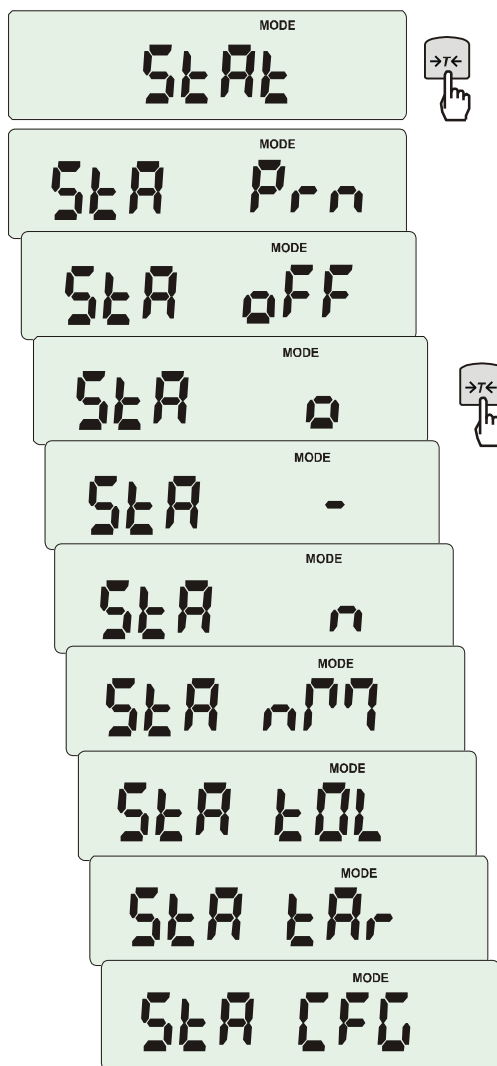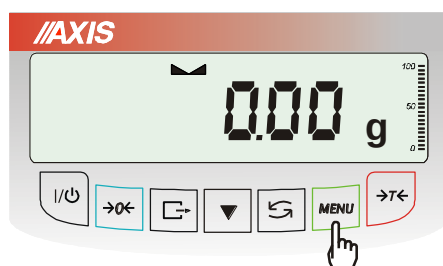

Press MENU key.

When StAt is displayed press →T← key.

The following options are displayed:

- StA Prn – monitoring and printout of statistical data,
- StA oFF – deactivate function,
- StA □ – activate function, work with printout of chosen weighting results,
- StA -- activate function, work without printout,
- StA n – maximal samples value,
- StA nM – inscribing nominal value for statistics,
- StA tOL – inscribing tolerance in %,
 - Auto – Automatic work (samples are confirmed after loading the scale and indication stabilization.),
 - ManuAL – manual work (confirmation is made by pressing →T← key).
- StA CFG – function configuration:
- out – exit from function.

Remember first to inscribe nominal weight value and tolerance (mentioned above).

After that, push →T← key when StA o is displayed.

Put on successive objects on the pan (remove after indication stabilization) in order to add them to measurements register.

In order to obtain printed statistical results from measurements series press MENU key and →T← key when StAt is displayed and later StA Prn.

After printout two options are enabled:

- rESet – erasing results,
- Contin – continuation.

Pressing  key printouts estimated values and histogram :

Nominal - nominal value,

Tolerance - accepted value in percentage.

N - number of sample

IN TOL. – number of samples in toleranc

-TOL – amount of measurements

under allowable lower value

+TOL – amount of measurements above

allowable upper value

TOTAL - sum of weights of all n samples

AVERAGE – average weight as (Total)/n

MIN – minimum weight in n samples

MAX– maximum weight in n samples

ST. DEV. – standard deviation

ST. DEV.% – standard deviation percentage

To finish work with this function and

zeroing result register press MENU

key and then when StAt. and Sta oFF is

displayed press →T← button.

Statistics function cooperation with computer and

printer: Scale can be equipped with two serial ports

marked as RS232C-I (computer) and RS232C-II

(printer). After each data printout by printer identical

set of data is sent to computer. After sending by

computer initialization signal S A CR LF

(53h 49h 0Dh 0Ah) the scale sends to computer

statistic data enclosed in histogram.

```

----- STATISTICS -----
NOMINAL :    13.20 g
TOLERANCE:    20 %
MAX. N :      500

```

NO.	SAMPLE	TOL-	NOM	TOL+
1	12.31 g	*		
2	12.78 g	*		
3	12.74 g	*		
4	12.75 g	*		
5	13.26 g		*	
6	13.22 g		*	
7	13.21 g		*	
8	13.22 g		*	
9	13.24 g		*	
10	13.70 g		*	
11	13.70 g		*	
12	13.72 g		*	
13	14.18 g		*	

```

N :      13
IN TOL. :    13
< TOL- :      0
> TOL+ :      0
TOTAL :    172.03 g
AVERAGE :    13.23 g
MAX :    14.18 g
MIN :    12.31 g
MAX-MIN :    1.87 g
ST.DEV. :    0.489 g
ST.DEV.% :    3.69 %

```

```

----- HISTOGRAM -----

```

```

<TOL-  0 |
        0 |
        0 |
        0 |
        1 |
        3 |
        5 |
        3 |
        1 |
        0 |
        0 |
        0 |
        0 |
>TOL+  0 |


```

```

SIGN. : .....

```

18.13 Function for summing recipe ingredients (*rECIPE*)

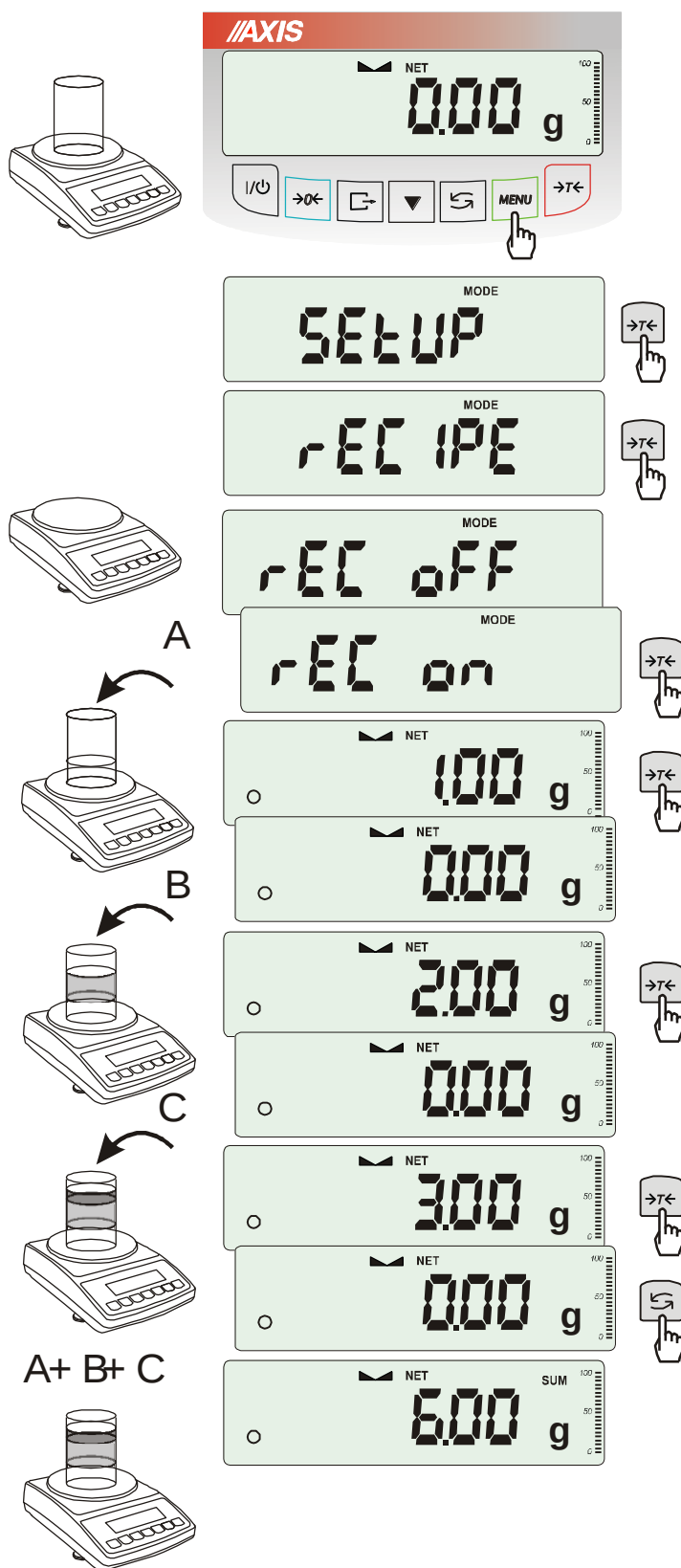

The function allows for separate weighing of several ingredients in one container with the possibility of reading current sum of all weighed ingredients.

The function includes the following options:

- rEC oFF* – leave the function with the possibility of read sum mass,
- rEC on* – start recipe weighing,
- rEC Con* – continue previous recipe,
- out* – exit without changes.

When proceeding with recipe, successive ingredients (A, B, C, etc.) are weighed each time starting from zero indication, which is obtained after scale taring.

If several ingredients are weighed, their sum mass can be read (despite several taring). For this purpose use **↺** or *rEC oFF* option.

Using **↺** once again enables fast return to recipe.



In order to turn off *rECIPE* function press **MENU** key and then using **→T←** key choose *rECIPE* and *rEC oFF*.

Comments:

o indicator on the left side of scale display shows *rECIPE* function is active.

SUM indicator shown when *rEC oFF* option is used, disappears after using **→T←** key.

18.14 Density determination (dEnSlty)

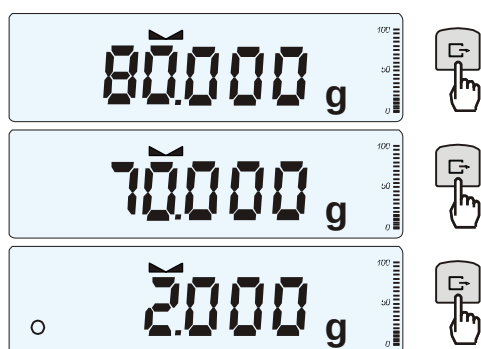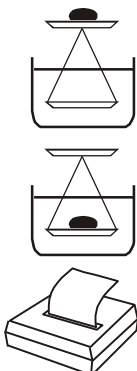

Solids density determination

This function calculates material density basing on its weight in air and in water using the formula below:

$$\rho = \frac{m_1}{m_1 - m_2} * \rho_L$$

where, m_1 – weigh in air
 m_2 – weight in water
 ρ_L – density of liquid

If distilled water (H_2O) or ethanol (*EthAnOL*) is used, enter its exact temperature (accurate to 0,5°C) – the balance will calculate its density automatically.

To enter the value use the following keys:

- ▼ - digit increase,
- ⇨ - decimal point,
- T← - next digit,
- MENU - end.

When using liquid other than distilled water or ethanol, choose *othEr* option and enter its density according to its temperature.


Phase I: measurement in air.
 Phase II: measurement in liquid.

To print measurement result and begin next measurement press ⇨ key.


To print a density determination report after all necessary measurements, connect a printer to the balance and press  key. A sample for solid mass density determination is shown below:

-----DENSITY-----	
Mass in the air	= ...
Mass in the liquid	= ... g
Mass density	= ... g/cm ³
Density with comp.	= ... g/cm ³
Water density	= ... g/cm ³
Water temperature	= ... °C

Operation sequence (weighing in air and in liquid) for below-balance weighing:

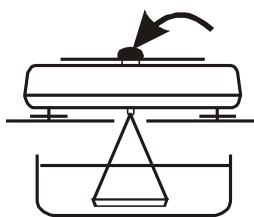

Phase I: measurement in air.

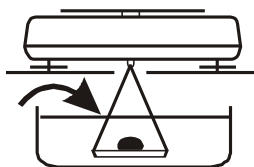

Phase II: measurement in liquid.

Liquid density determination



This function determines liquid density basing on plunger weight in air and in examined liquid with known volume, using the formula below:

$$\rho = \frac{m_1 - m_2}{V}$$

where

m_1 – plunger weight in air

m_2 – plunger weigh in a liquid

V – plunger volume

Plunger volume is stored on its hanger.

To enter the value use the following keys:

▼ - digit increase,

↔ - decimal point,

→T← - next digit,


MENU - end.

Phase I: measurement in air.

Phase II: measurement in liquid.

To print measurement result and begin next measurement press ↔ key.

18.15 Paper weight determination function (PAPeR) - optional

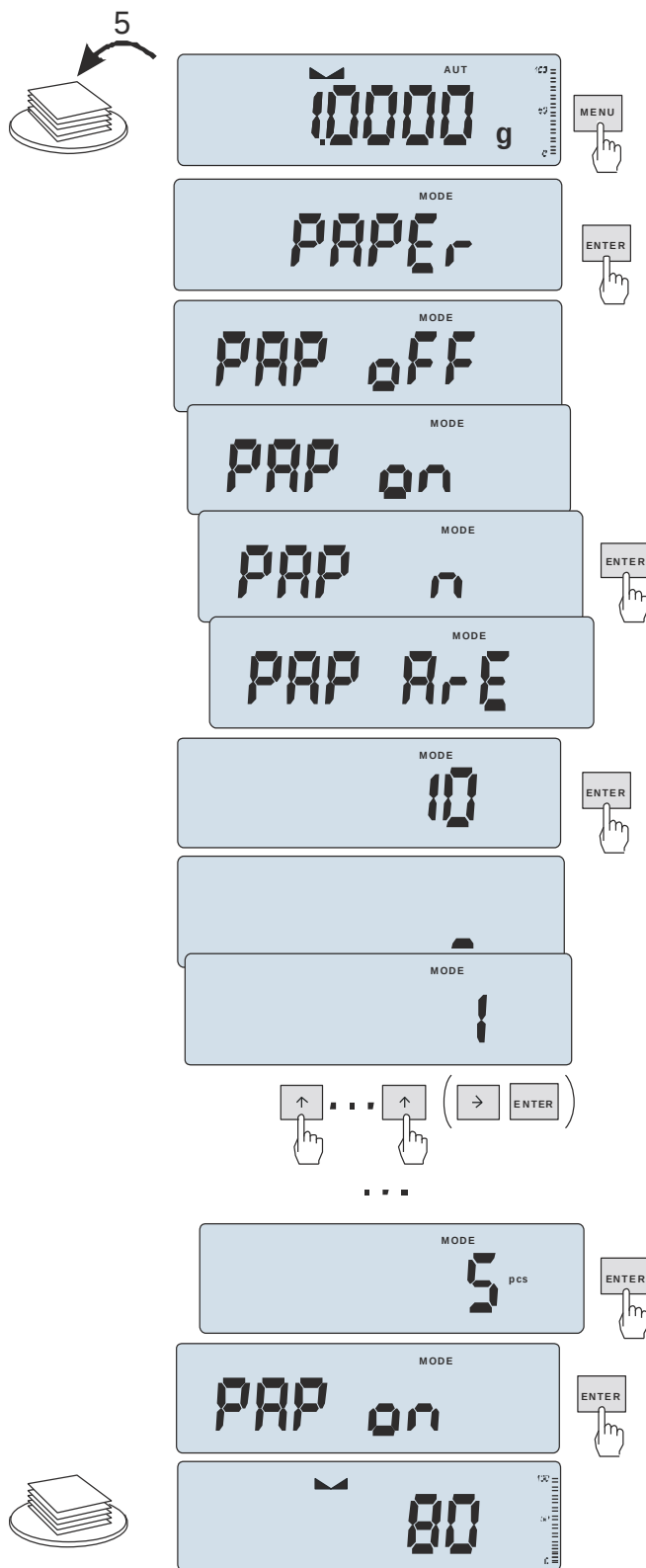

This function allows you to calculate the weight of 1 m² of paper based on a sample of several sections with a known area.

This function allows you to calculate the mass of 1 m² of paper based on a sample of several paper cutouts of known area. Tare the scale using the $\rightarrow T \leftarrow$ key.

Load a sample of one or more paper cutouts onto the weighing pan (ensure that the total load is no less than 100 scale reading units).

Press the MENU key to access the function menu. Select the PAPER function.

The display will show the following:

- *PAP OFF* – disables the function
- *PAP ON* – measures the basis weight in g/m²
- *PAP n* – enters the number of sections placed on the pan
- *PAP ArE* – enters the area of a single section in m²

To enter *PAP n* and *PAP ArE*, use the following keys::

$\rightarrow 0 \leftarrow$ - increase digit,

$\rightarrow T \leftarrow$ - move to the next digit,

MENU – end entry.

Press the $\rightarrow T \leftarrow$ key when the display shows *PAP ON*.

The scale will indicate the paper weight in g/m², indicated by the "o" symbol on the left side of the display.

To exit the function, press the MENU key, and then use the $\rightarrow T \leftarrow$ key to select *PAPER* and *PAP OFF*.

Uwaga:

1. Komunikat "PAP Err" oznacza, że wpisano niewłaściwe wartości *PAP n* lub *PAP ArE*.

19. Troubleshooting and maintenance

1. The scale should be kept clean.
2. Take care that no dirt is between pan and casing of the scale. If dirt is noticed take the pan off (lift it up). Clean dirt and then put the pan on.
3. In case of improper operation caused by a short-lasting lack of power supply, switch the scale off by unplugging it from the mains, and then after several seconds switch it on.
4. Every repairs performed by unauthorized persons are forbidden
5. To repair the scale, please contact nearest service centre. The list of authorised service centres is given in guarantee card and on www.axis.pl website.
6. Scales can be sent for repair as messenger delivery only in original package. For transportation scale pan have to be protected against accidental pressing. If not, there is a risk of damaging the scale and loosing guarantee.

Failure messages:

Message	Possible cause	Recommendation
C-1 ... 6 (more than 1min.)	negative result in one of autotests	contact service centre if the message remains
scale is not weighing	protecting screw remains in the scale	remove protecting screw
L	no pan on the scale	put the pan on
	mechanical damage of scale sensor	contact service centre
H	overweight of the scale	take a load off the pan
	mechanical damage of the scale	contact service centre
UnLOAD	load left on the pan during start-up	take a load off the pan
SERVICE	mechanical damage	contact service centre
 indicator does not work	unstable scale position, ground vibration, air flows	locate the scale in place where stable results are maintained
	damage of the scale	contact service centre
-----	taring not finished	contact service centre
Err t	przekroczenie dopuszczalnego zakresu temperatur pracy wagi	przenieść wagę do miejsca o właściwej temperaturze
Err C09	zakłócony przebieg lub wadliwy wynik kalibracji wewnętrznej	powtórzyć kalibrację wewnętrzną, jeśli komunikat pozostaje, zgłosić wagę do serwisu

Notes